GDDAQ

发布 *beta*

Hongyi Wu(吴鸿毅)

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content:

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Event header as the first 4 words
RevD(12-bit,100MHz),RevF(14-bit,100MHz)

Index	Data					
	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
0	Finish Code	Event Length	Header Length	CrateID	SlotID	Chan#
1	[31:0]					
	EVTTIME_LO[31:0]					
2	[31]	[30:16]			[15:0]	
	CFD forced trigger bit	CFD Fractional Time[14:0] x 32768			EVTTIME_HI[15:0]	
3	[31]	[30:16]			[15:0]	
	Trace Out-of-Range Flag	Trace Length			Event Energy	

(EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²+(CFD_Fractional_Time[14:0]/32768))x10ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

Event header as the first 4 words RevF(12/14/16-bit,250MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31]	[30]	[29:16]	[15:0]		
	CFD forced trigger bit	CFD trigger source bit	CFDFractionalTime[13:0]x16384	EVTTIME_HI[15:0]		
3	[31]	[30:16]	[15:0]			
	Trace Out-of-Range Flag	Trace Length	Event Energy			

((EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²)x2 - CFD trigger source bit + (CFD_Fractional_Time[13:0]/16384))x4ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

CFD trigger source bit:

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

Event header as the first 4 words RevF(12/14-bit,500MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31:29]	[28:16]	[15:0]			
	CFD trigger source bits[2:0]	CFD Fractional Time[12:0] x 8192	EVTTIME_HI[15:0]			
3	[31]	[30:16]	[15:0]			
	Trace Out-of-Range Flag	Trace Length	Event Energy			

(EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²)x10 + ((CFD_Fractional_Time[12:0]/8192)+ CFD trigger source bits[2:0]-1)x2ns

Finish Code: 0-good event,1- pileup event

CFD trigger source bits:

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

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
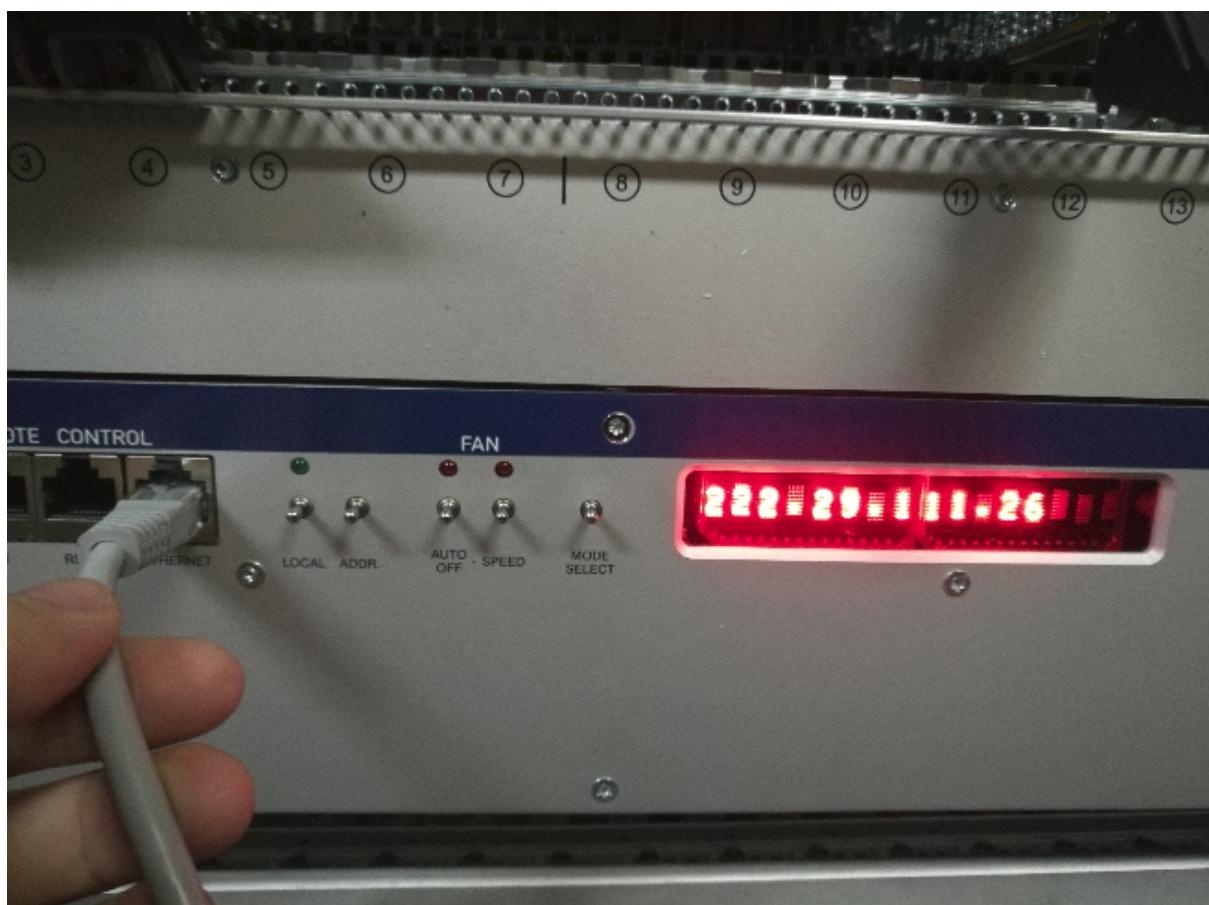
External Temperature Sensors and Voltage Inputs

1	2	3	4	5	6	7	8
66°F		66°F		66°F			

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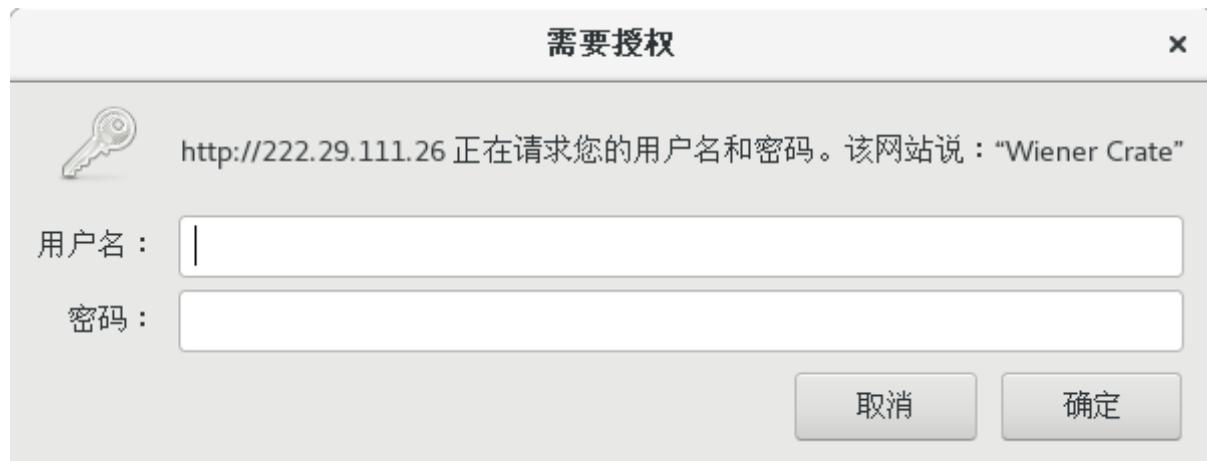The screenshot shows a web-based monitoring interface for the UEP6000/PL500 system. At the top, there's a header bar with the title "UEP6000/PL500 - Mozilla Firefox". Below it is a navigation bar with a back button, forward button, address bar containing "222.29.111.26", a search bar with placeholder "搜索", and various browser control icons.

Below the navigation bar, there are four buttons: "MAIN POWER", "VME SYSRESET", "FAN SLOWER", and "FAN FASTER". To the right of these buttons, the text "W-IE-NE-R" is displayed.

The main content area is divided into sections:

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
7.3.4 CFD

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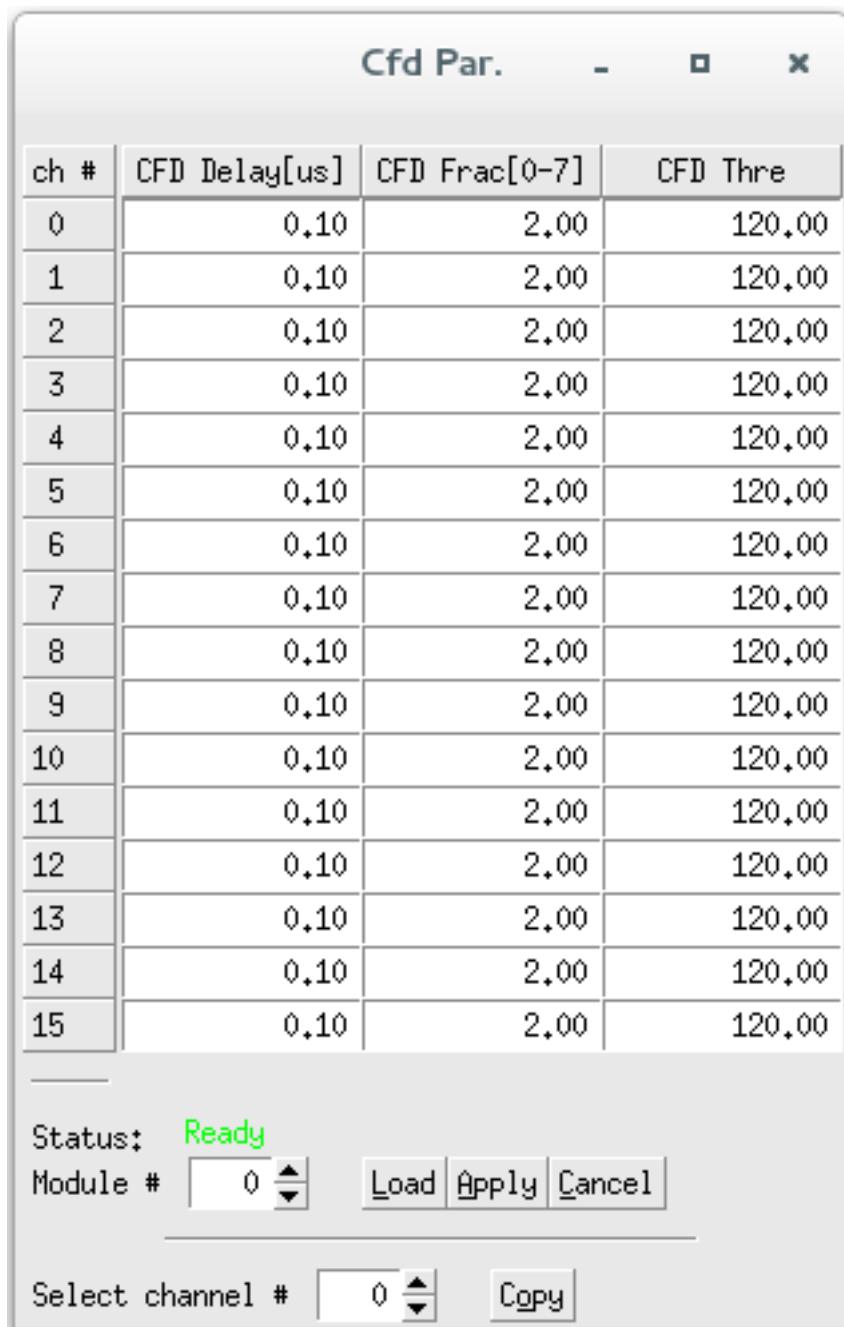TODO

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The following CFD algorithm is implemented in the signal processing FPGA of the 100 MHz(Rev. B, C, D and F) and 250 MHz(Rev. F) Pixie-16 modules.

Assume the digitized waveform can be represented by data series $\text{Trace}[i]$, $i = 0, 1, 2, \dots$. First the fast filter response(FF) of the digitized waveform is computed as follows:

$$\text{FF}[i] = \sum_{j=i-(FL-1)}^i \text{Trace}[j] - \sum_{j=i-(2 \times FL+FG-1)}^{i-(FL+FG)} \text{Trace}[j]$$

Where FL is called the fast length and FG is called the fast gap of the digital trapezoidal filter. Then the CFD is computed as follows:

$$CFD[i + D] = FF[i + D] \times (1 - w/8) - FF[i]$$

Where D is called the CFD delay length and w is called the CFD scaling factor($w=0, 1, \dots, 7$).


The CFD zero crossing point(ZCP) is then determined when $CFD[i] \leq 0$ and $CFD[i + 1] < 0$. The timestamp is latched at Trace point i , and the fraction time f is given by the ratio of the two CFD response amplitudes right before and after the ZCP.

$$f = \frac{CFDout1}{CFDout1 - CFDout2}$$

Where CFDout1 is the CFD response amplitude right before the ZCP, and CFDout2 is the CFD response amplitude right after the ZCP(subtraction is used in the denominator since CFDout2 is negative). The Pixie-16 DSP computes the CFD final value as follows and stores it in the output data stream for online or offline analysis.

$$CFD = \frac{CFDout1}{CFDout1 - CFDout2} \times N$$

Where N is scaling factor, which equals to 32768 for 100 MHz modules and 16384 for 250 MHz modules, respectively.

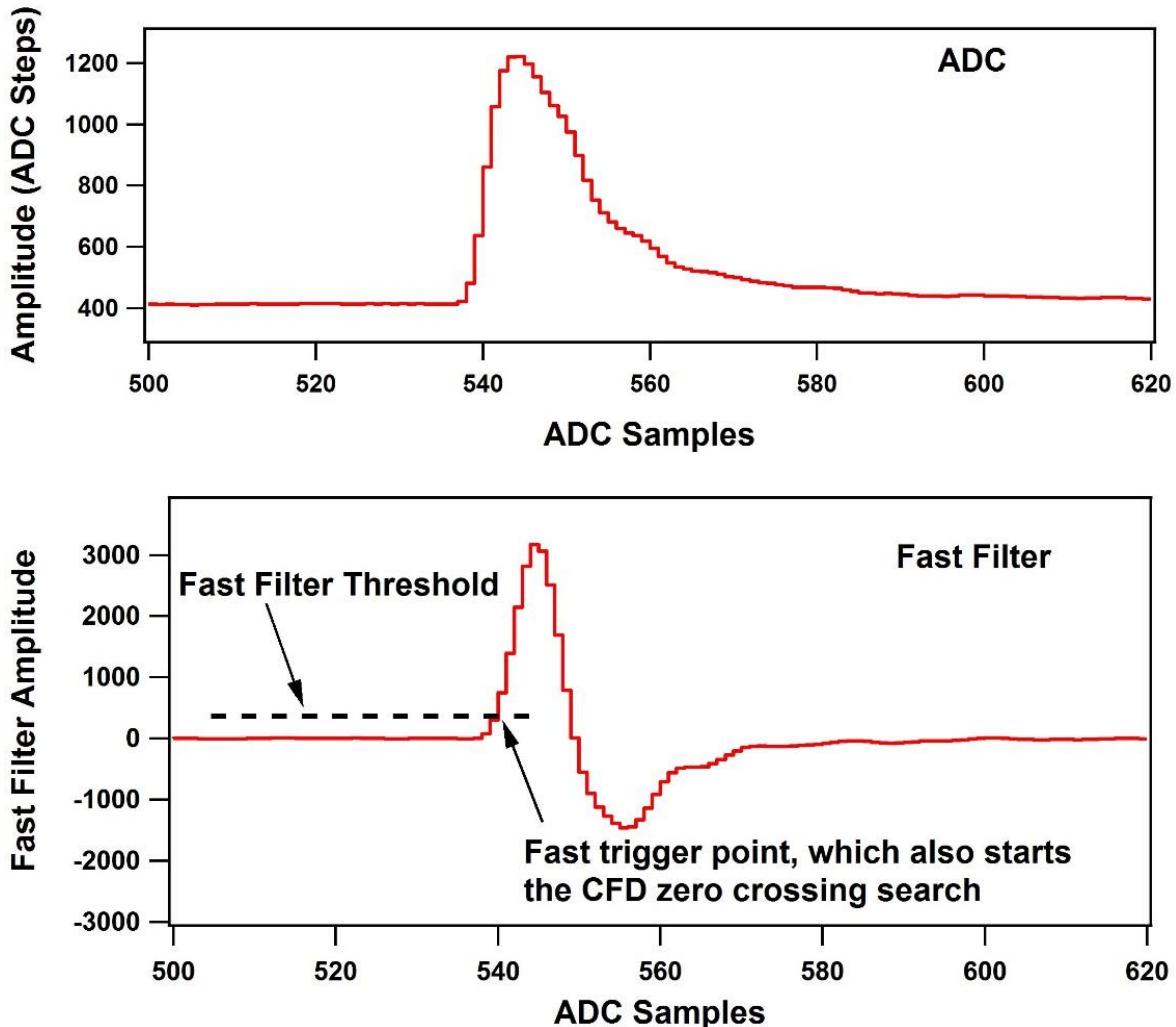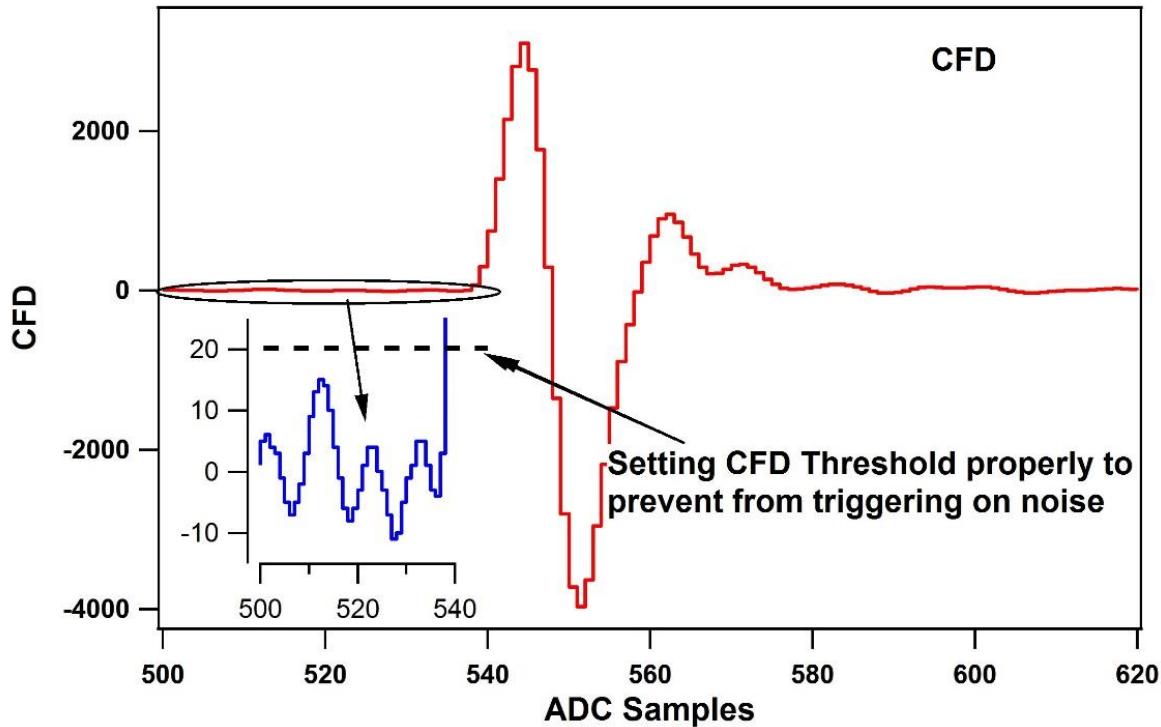
Figure shows a sample ADC trace, its fast filter response and its CFD response, respectively.

The top figure shows a raw ADC trace. After computing the fast filter response on the raw ADC trace using Equation $FF[i]$, the fast filter response is compared against the fast filter threshold as shown in the middle figure. The ADC sample where the fast filter response crosses the fast filter threshold is called the fast trigger point, which also starts the search for the CFD zero crossing point.

The CFD response is computed using Equation $CFD[i + D]$ and is shown in the bottom figure(for actual implementation in the firmware, the fast filter response FF is delayed slightly before being used for computing the CFD response so that there are sufficient number of CFD response points to look for the zero crossing point after the fast trigger). To prevent premature CFD trigger as a result of the noise in the CFD response before the actual trigger, a DSP parameter called CFDThresh is used to suppress those noise-caused zero crossing. However, if a zero crossing point cannot be found within a certain period after the fast trigger (typically 32 clock cycles), e.g., due to unnecessarily high CFDThresh, a forced CFD Trigger will be issued and a flag will be set in an event header word to indicate that the recorded CFD time for this event is invalid.

However, the event will still have a valid timestamp which is latched by the fast filter trigger when fast filter crosses over the trigger threshold. The aforementioned CFD parameters correspond to the following DSP parameters.

CFD Parameters	DSP Parameters
FL	FastLength
FG	FastGap
Fast Filter Threshold	FastThresh
D	CFDDelay
W	CFDScale (valid values: 0, 1, 2, ... and 7)
CFD Threshold	CFDThresh

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In the 250 MHz Pixie-16 modules, the event timestamp is counted with 125 MHz clock ticks, i.e., 8 ns intervals, and two consecutive 250 MHz ADC samples are captured in one 8 ns interval as well.

The CFD trigger also runs at 125 MHz, but the CFD zero crossing point is still reported as a fractional time between two neighboring 250 MHz ADC samples, which are processed by the FPGA in one 125 MHz clock cycle.

However, the CFD zero crossing point could be in either the odd or even clock cycle of the captured 250 MHz ADC waveforms.

Therefore, the firmware outputs a “CFD trigger source” bit in the output data stream to indicate whether the CFD zero crossing point is in the odd or even clock cycle of the captured 250 MHz ADC waveforms.

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In the 100 MHz Pixie-16 modules, event timestamp, CFD trigger, and ADC waveform capture are all carried out with the same 100 MHz clock. So there is no need to report “CFD trigger source” for the 100 MHz Pixie-16 modules.

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The CFD algorithm discussed in the previous section for the 100 MHz and 250 MHz Pixie-16 modules can also be written in the following format:

$$CFD(k) = w \cdot \left(\sum_{i=k}^{k+L} a(i) - \sum_{i=k-B}^{k-B+L} a(i) \right) - \left(\sum_{i=k-D}^{k-D+L} a(i) - \sum_{i=k-D-B}^{k-D-B+L} a(i) \right)$$

Where $a(i)$ is the ADC trace data, k is the index, and w , B , D , and L are CFD parameters.

The CFD algorithm implemented in the 500 MHz Pixie-16 modules is special when compared to the one implemented in the 100 MHz and 250 MHz Pixie-16 modules in terms of the ability to adjust parameters w , B , D , and L .

The reason for this is that in the 500 MHz Pixie-16 modules, ADC data that come into the FPGA at the speed of 500 MHz is first slowed down with a ratio of 1:5, in other words, the FPGA captures 5 ADC samples at the rate of 100 MHz, i.e., every 10 ns. The FPGA then tries to find the CFD trigger point between any two adjacent 2-ns ADC samples within that 10 ns by first building sums of ADC samples and then calculating differences between delayed and non-delayed sums until the zero crossing point is found. However, in the 500 MHz Pixie-16 modules, the FPGA does not have enough resources to build sums for 5 ADC samples in parallel with variable delays. Therefore, the CFD algorithm for the 500 MHz modules was implemented using a set of fixed CFD parameters as shown in Table *Fixed CFD Parameter Values for 500 MHz Pixie-16 Modules*. Tests show these fixed parameters give best performance for LaBr₃(Ce) detectors.

CFD Parameters	Fixed Values for 500 MHz Modules
w	1
B	5
D	5
L	1

The CFD time given by the 500 MHz Pixie-16 modules consists of two parts: a shift within the 5 ADC samples and a fractional time between two ADC samples where the CFD zero crossing occurred. The shift within the 5 ADC samples is reported as the 3-bit CFD trigger source[2:0] is defined as follows.


CFD Trigger Source [2:0]	Zero Crossing Point (ZCP) Location
000	ZCP occurred between the 5th ADC sample of the previous 5-sample group and the 1st ADC sample of the current 5-sample group
001	ZCP occurred between the 1th ADC sample of the current 5-sample group and the 2nd ADC sample of the current 5-sample group
010	ZCP occurred between the 2nd ADC sample of the current 5-sample group and the 3rd ADC sample of the current 5-sample group
011	ZCP occurred between the 3rd ADC sample of the current 5-sample group and the 4th ADC sample of the current 5-sample group
100	ZCP occurred between the 4th ADC sample of the current 5-sample group and the 5th ADC sample of the current 5-sample group
101	Not used
110	Not used
111	CFD trigger is forced, so CFD time is invalid

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$$CFD = \frac{CFDout1}{CFDout1 - CFDout2} \times 8192$$

7.3.5 QDC

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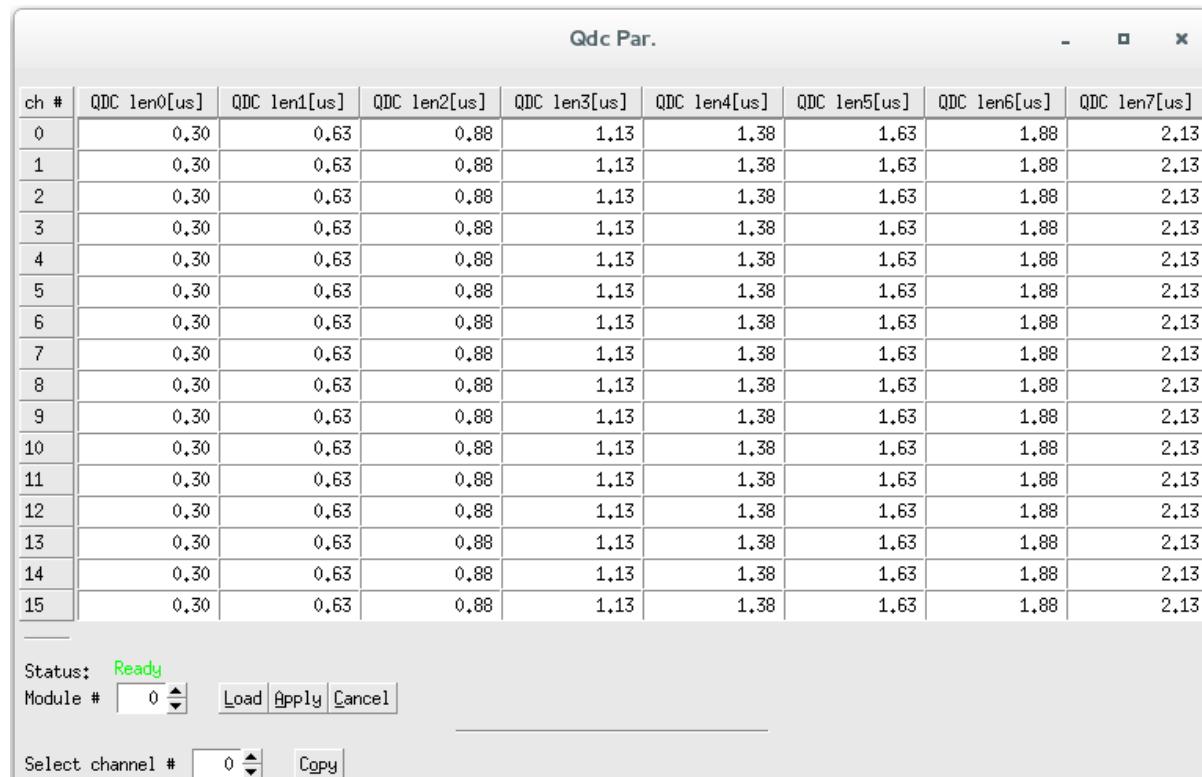
Eight QDC sums, each of which can have different lengths, are computed in the Signal Processing FPGA of a Pixie-16 module for each channel and the sums are written to the list mode output data stream if the user requests so.


The recording of QDC sums starts at the waveform point which is *Pre-trigger Trace Length* or *Trace Delay* earlier than the trigger point, which is either the CFD trigger or channel fast trigger depending on whether or not CFD trigger mode is enabled.

The eight QDC sums are computed one by one continuously, but they are not overlapping. The recording of QDC sums ends when the eight intervals have all passed.

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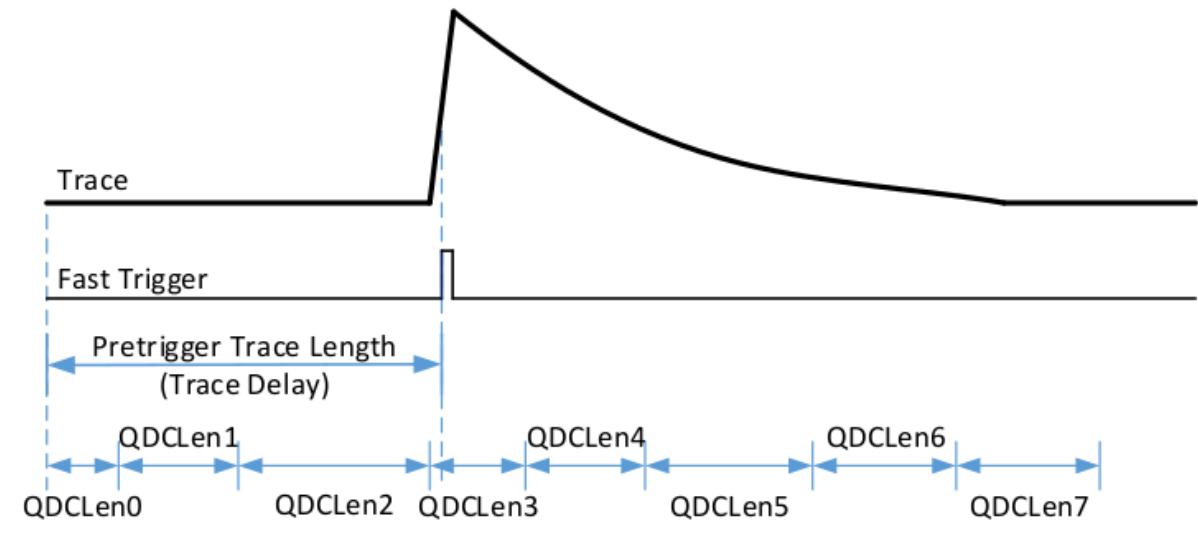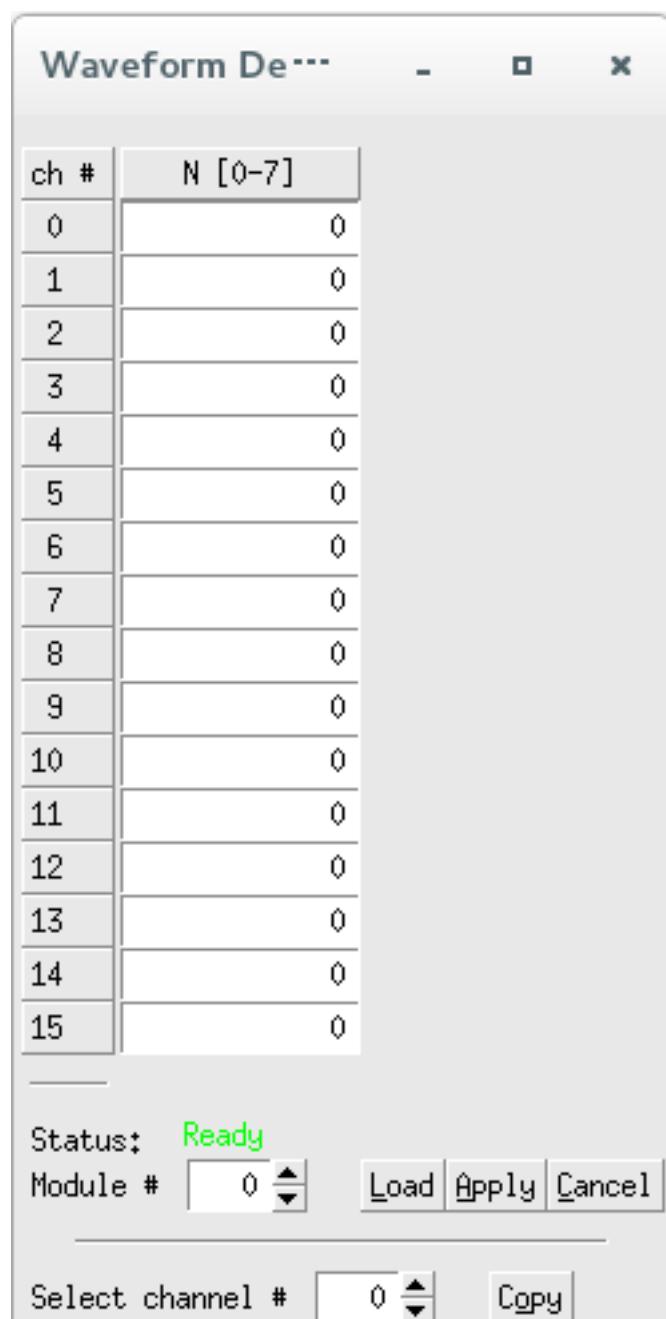
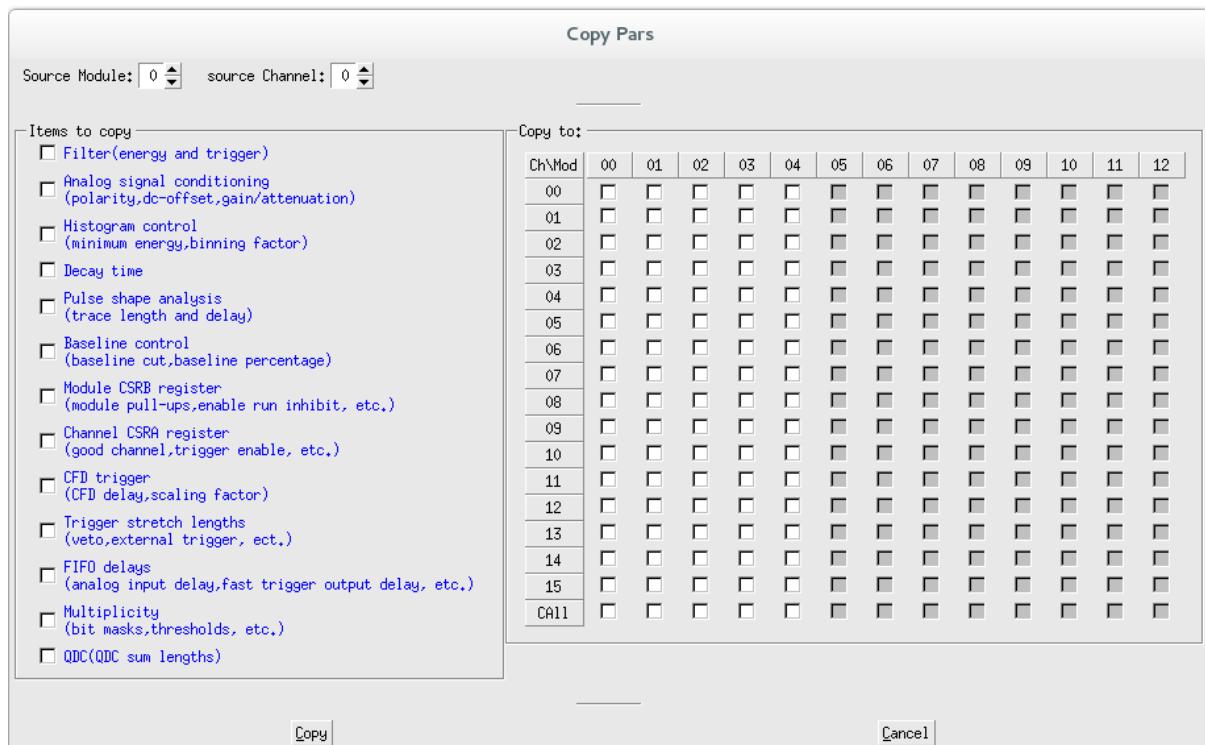
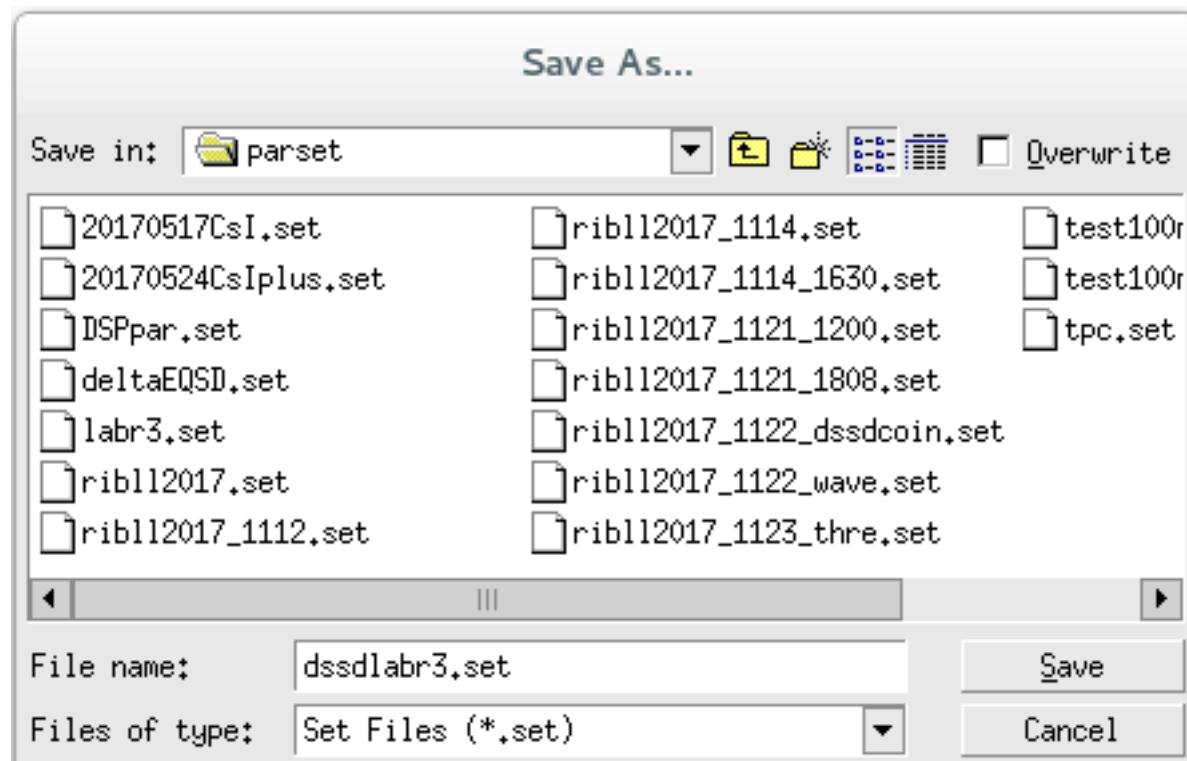
7.4.1 Module Variables

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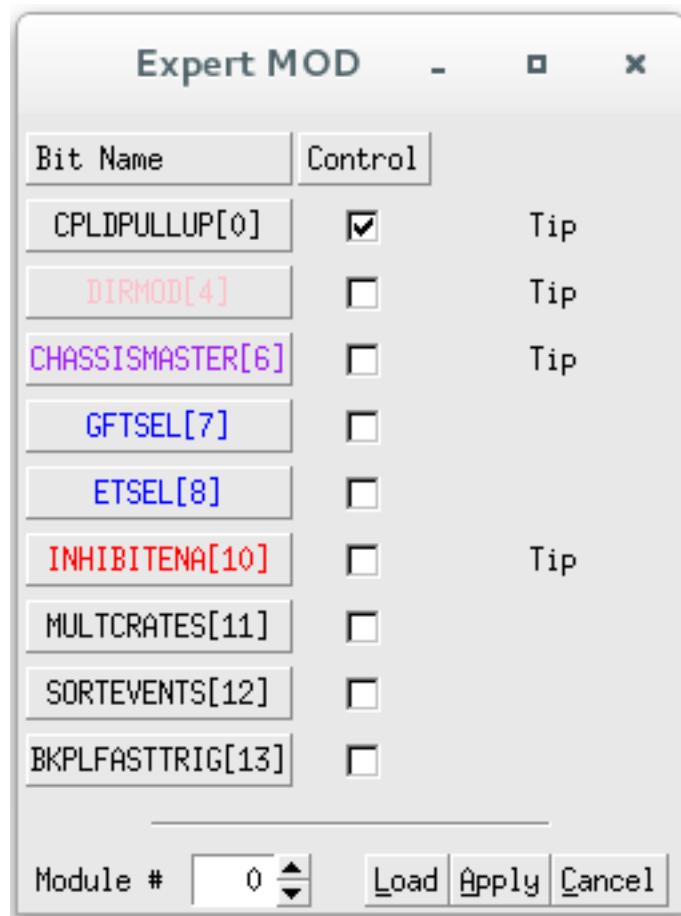In addition to distributing the global clock signal, the Pixie-16 rear I/O trigger module can also share global triggers and run synchronization signals. The global trigger signals include the global validation trigger and global fast trigger, plus the Pixie-16 FPGA data storage buffers' full flag signal. The run synchronization signals include synchronous run start and stop signals that can be shared among multiple crates.

In order to enable the distribution of such global triggers and run synchronization signals, certain Pixie-16 parameters have to be set properly. The parameter that controls the trigger distribution and run synchronization is the Module Control Register B (ModCSRB).

ModCSRB is a 32-bit parameter with each of 32 bits controlling different operation modes of the Pixie-16 module.

注解：Trigger Distribution and Run Synchronization

For the System Director module that is installed in the Master crate, bits 0, 4, 6 and 11 of ModCSRB should be set to 1 (checked & enabled).

For the Crate Master module that is installed in the Slave crate, bits 0, 6 and 11 of ModCSRB should be set to 1 (checked & enabled).

For the General modules that are installed in both the Slave crate and Master crate, bit 11 of ModCSRB should be set to 1 (checked & enabled).

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Module Control Register B affecting the module as a whole.

- bit 0 - MODCSRB_CPLDPULLUP

- Enable pullups for PXI trigger lines on the backplane through an onboard CPLD.
 - With the pullups, those PXI trigger lines default to logic high state.
 - Only when one module actively pulls a line to logic low state will such a line be in the low state.
 - Therefore signals transmitted over those PXI trigger lines are actively low signals.
 - **Note: enable this bit only for one module per crate (e.g. the crate master module)**
 - **bit 4 - MODCSRB_DIRMOD**
 - Set this module as the Director module so that it can send triggers, trace and header DPM full signal and run synchronization signal to all crates through the rear I/O trigger modules.
 - Here triggers include fast trigger and validation trigger
 - **Note: enable this bit only for one module among all crates (e.g. the system director module in multi-crate configuration)**
 - **bit 6 - MODCSRB_CHASSISMASTER**
 - Set this module as the chassis master module so that it can send triggers, trace and header DPM full signal and run synchronization signal to the backplane of the local crate.
 - Here triggers include fast trigger and validation trigger
 - **Note: enable this bit only for one module per crate(e.g. the crate master module)**
 - **bit 7 - MODCSRB_GFTSEL**
 - Select external fast trigger source(=1: external validation trigger, =0: external fast trigger, in case these two signals are swapped at the Pixie-16 front panel input connectors)
 - **bit 8 - MODCSRB_ETSEL**
 - Select external validation trigger source(=1: external fast trigger, =0: external validation trigger, in case these two signals are swapped at the Pixie-16 front panel input connectors)
 - **bit 10 - MODCSRB_INHIBITENA**
 - Enable(=1) or disable(=0) the use of external INHIBIT signal.
 - When enabled, the external INHIBIT signal in the logic high state will prevent the run from starting until this external INHIBIT signal goes to logic low state.
 - **bit 11 - MODCSRB_MULTCRATES**
 - Set this module to run in the multi-crate mode(=1) or in the local-crate mode(=0).
 - If the module is running in multi-crate mode, it will use the trace and header DPM full signal and run synchronization signal that are generated and distributed among multiple crates.
 - If the module is running in local-crate mode, it will use the trace and header DPM full signal and run synchronization signal generated in the local crate.
 - **bit 12 - MODCSRB_SORTEVENTS**
 - Sort(=1) or do not sort(=0) events from all 16 channels of a Pixie-16 module based on the timestamps of the events, before storing the events in the external FIFO.
 - Note: all 16 channels must have the same DAQ parameters setting to use this feature
 - **bit 13 - MODCSRB_BKPLFASTTRIG**
 - Enable(=1) or disable(=0) the sending of 16 local fast triggers to the 16 lines on the backplane of the crate.
 - **Note: only one module can enable this option in each PCI bus segment of a crate(not limited to the crate master module, e.g. any module in each PCI bus segment)**
-

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Channel Control Register A affecting each channel individually


- **bit 0 - CCSRA_FTRIGSEL**
 - Channel fast trigger selection(=1: module fast trigger from the System FPGA; =0: the selection depends on the value of another bit CCSRA_GROUPTRIGSEL – if CCSRA_GROUPTRIGSEL = 1, select the channel validation trigger from the System FPGA, and if CCSRA_GROUPTRIGSEL = 0, select this channel’s local fast trigger)
- **bit 1 - CCSRA_EXTTRIGSEL**
 - Module validation trigger selection(=1: module gate input from the Pixie-16 front panel Module Gate LVDS connector; =0: module validation trigger from the System FPGA)
- **bit 2 - CCSRA_GOOD**
 - Set this channel as a Good channel(=1) or a not Good channel(=0).
 - **When a channel is set to be a not Good channel, it still generates local fast triggers, which could be used in multiplicity computation, etc., but this channel will not record list mode data or MCA data, and will not update its baseline value**
- **bit 3 - CCSRA_CHANTRIGSEL**
 - Channel validation signal selection(=1: channel gate input from the Pixie-16 front panel Channel Gate LVDS connector; =0:channel validation trigger from the System FPGA)
- **bit 4 - CCSRA_SYNCDATAACQ**
 - Choose the level of synchronous data acquisition for this channel(=1: stops taking data when the trace or header DPM for any channel of any Pixie-16 module in the system is full; =0: stops taking data only when the trace or header DPM for this channel of this Pixie-16 module is full)
- **bit 5 - CCSRA_POLARITY**
 - Choose this channel’s input signal polarity(=1: invert input signal’s polarity; =0: do not invert input signal’s polarity).
 - **Please note in Pixie-16, signal processing requires positive rising input signal. So if input signal has a negative falling edge, it should be inverted by setting this CCSRA_POLARITY bit to 1**
- **bit 6 - CCSRA_VETOENA**
 - Enable(=1) or disable(=0) this channel’s veto.
 - If veto is enabled, this channel’s fast trigger will be vetoed by either the module veto signal(see bit 20 CCSRA_MODVETOSEL below) or channel veto signal(see bit 19 CCSRA_CHANVETOSEL below).
 - But if veto is disabled, this channel’s fast trigger will not be vetoed by either veto signal, even if either veto signal is present
- **bit 7 - CCSRA_HISTOE**
 - Enable(=1) or disable(=0) the histogramming of pulse energy values in the onboard MCA memory.
 - However, the current Pixie-16 firmware always histograms pulse energy values in the onboard MCA memory.
 - So this CCSRA_HISTOE is essentially not in use at the moment
- **bit 8 - CCSRA_TRACEENA**
 - Enable(=1) or disable(=0) trace capture in the list mode run for this channel
- **bit 9 - CCSRA_QDCENA**
 - Enable(=1) or disable(=0) QDC sums recording in the list mode run for this channel.

- There are a total of 8 QDC sums for each event
- **bit 10 - CCSRA_CFDMODE**
 - Enable(=1) or disable(=0) CFD trigger in the list mode run for this channel.
 - CFD trigger is used to latch sub-sample timing for the event time of arrival or timestamp
- **bit 11 - CCSRA_GLOBTRIG**
 - Enable(=1) or disable(=0) the requirement of module validation trigger for this channel.
 - If enabled, only when module validation trigger overlaps the channel fast trigger will the events be recorded for this channel
- **bit 12 - CCSRA_ESUMSENA**
 - Enable(=1) or disable(=0) the recording of raw energy sums and baseline values in the list mode run for this channel.
 - There are a total of three raw energy sums and one baseline value for each event.
 - **Please note the baseline value is stored in the format of 32-bit IEEE float point(IEEE 754)**
- **bit 13 - CCSRA_CHANTRIG**
 - Enable(=1) or disable(=0) the requirement of channel validation trigger for this channel.
 - If enabled, only when channel validation trigger overlaps the channel fast trigger will the events be recorded for this channel
- **bit 14 - CCSRA_ENARELAY**
 - Switch between two attenuations or gains for the input signal in this channel through an input relay(=1: close the input relay resulting in no input signal attenuation; =0: open the input relay resulting in a 1/4 input signal attenuation)
- **bit 15/16 - CCSRA_PILEUPCTRL/CCSRA_INVERSEPILEUP**
 - Control normal pileup rejection(bit 15) and inverse pileup rejection(bit 16) for list mode runs:
 - 1) Bits [16:15] = 00, record all events
 - 2) Bits [16:15] = 01, only record single events, i.e., reject piled up events
 - 3) Bits [16:15] = 10, record everything for piled up events, but will not record trace for single events even if trace recording is enabled, i.e., only record event header
 - 4) Bits [16:15] = 11, only record piled up events, i.e., reject single events
 - In all cases, if the event is piled up, no energy will be computed for such event
- **bit 17 - CCSRA_ENAENERGYCUT**
 - Enable(=1) or disable(=0) the “no traces for large pulses” feature.
 - If enabled, trace will not be recorded if the event energy is larger than the value set in DSP parameter EnergyLow
- **bit 18 - CCSRA_GROUPTRIGSEL**
 - Select channel fast trigger – this bit works together with the CCSRA_FTRIGSEL bit(bit 0): if CCSRA_FTRIGSEL=1, this CCSRA_GROUPTRIGSEL bit has no effect; if CCSRA_FTRIGSEL=0, then if CCSRA_GROUPTRIGSEL=1, select the channel validation trigger from the System FPGA, and if CCSRA_GROUPTRIGSEL=0, select this channel’s local fast trigger
- **bit 19 - CCSRA_CHANVETOSEL**
 - Channel veto signal selection(=1: channel validation trigger from the System FPGA; =0: channel gate input from the Pixie-16 front panel Channel Gate LVDS connector)
- **bit 20 - CCSRA_MODVETOSEL**

- Module veto signal selection(=1: module validation trigger from the System FPGA; =0: module gate input from the Pixie-16 front panel Module Gate LVDS connector)
- **bit 21 - CCSRA_EXTTSENA**
 - Enable(=1) or disable(=0) the recording of the 48-bit external clock timestamp in the event header during list mode run for this channel

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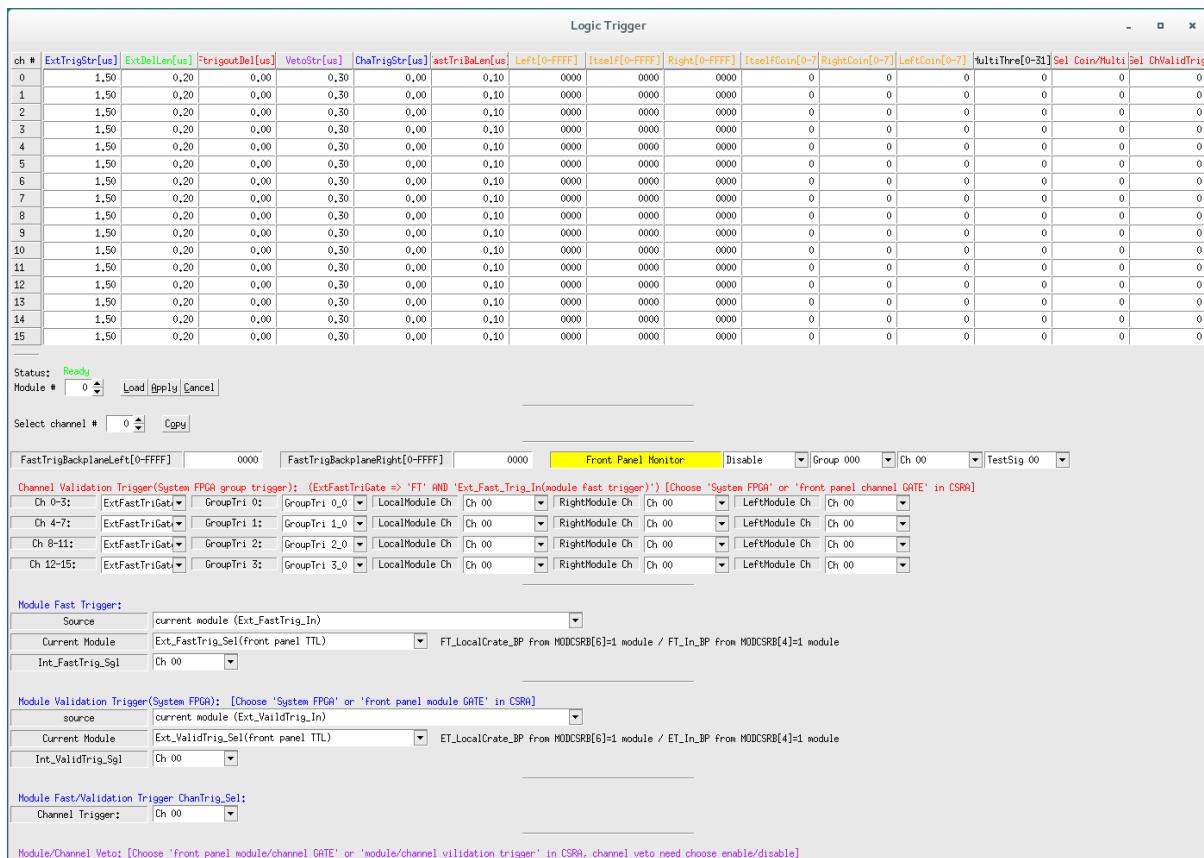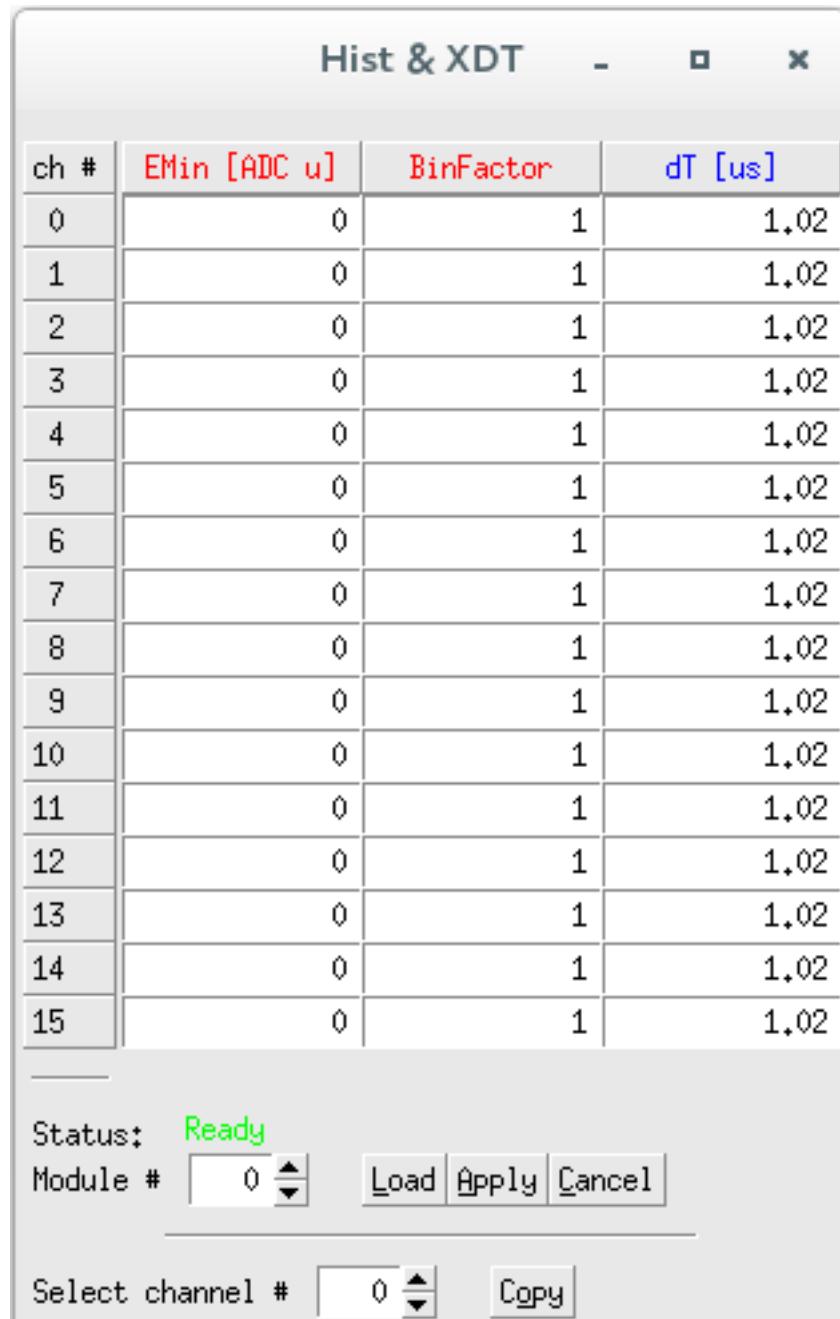
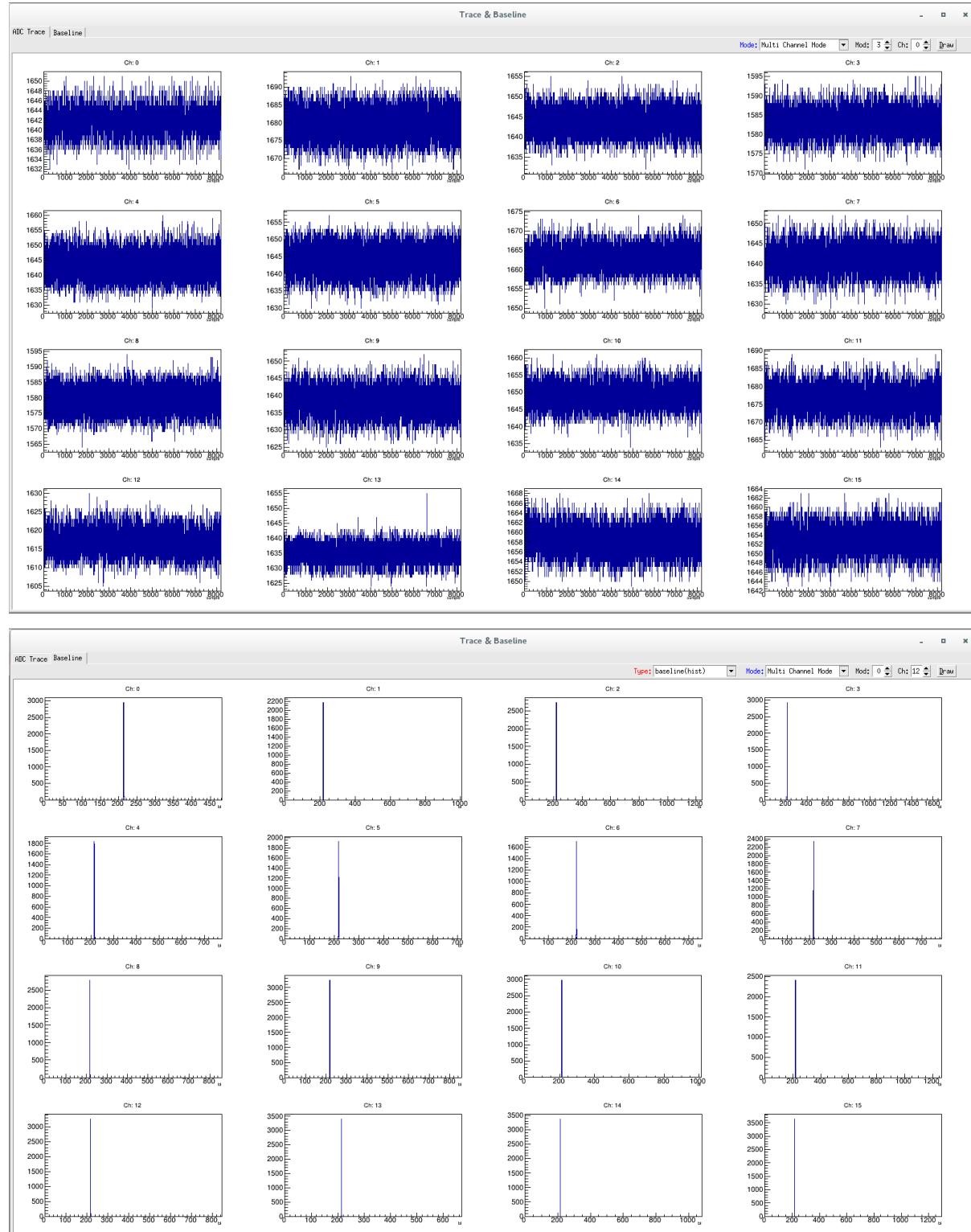
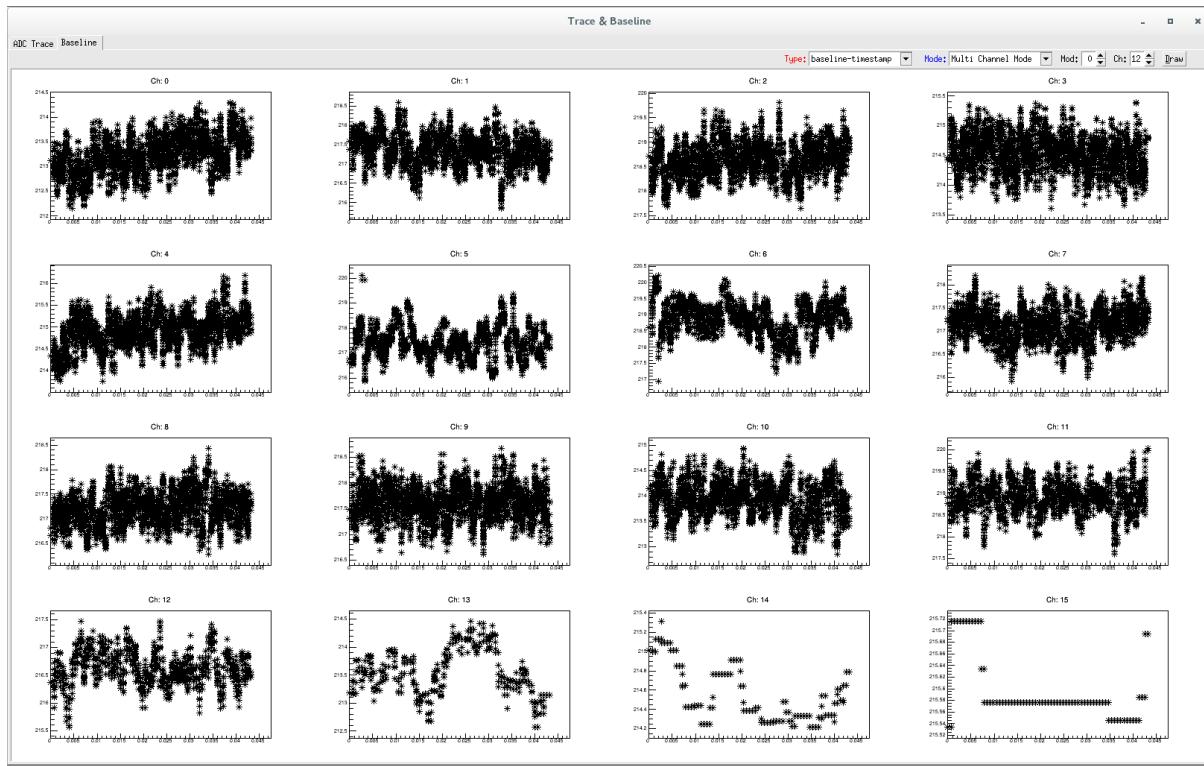
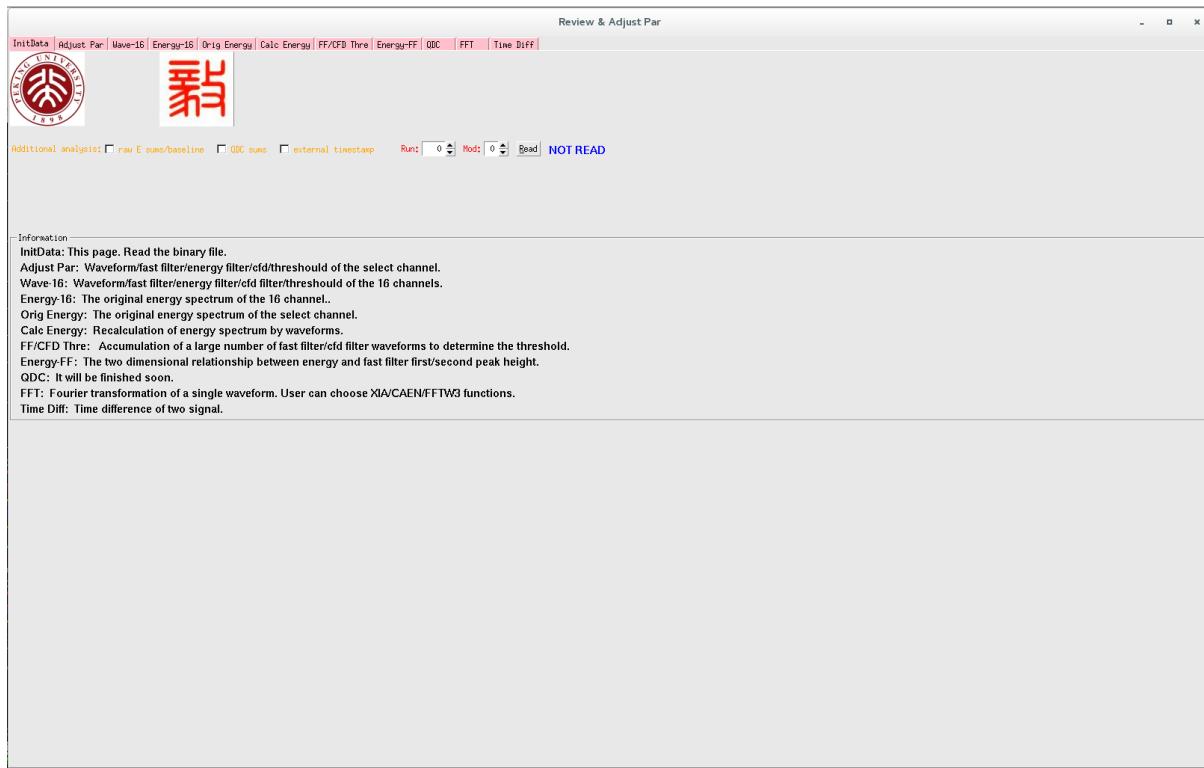
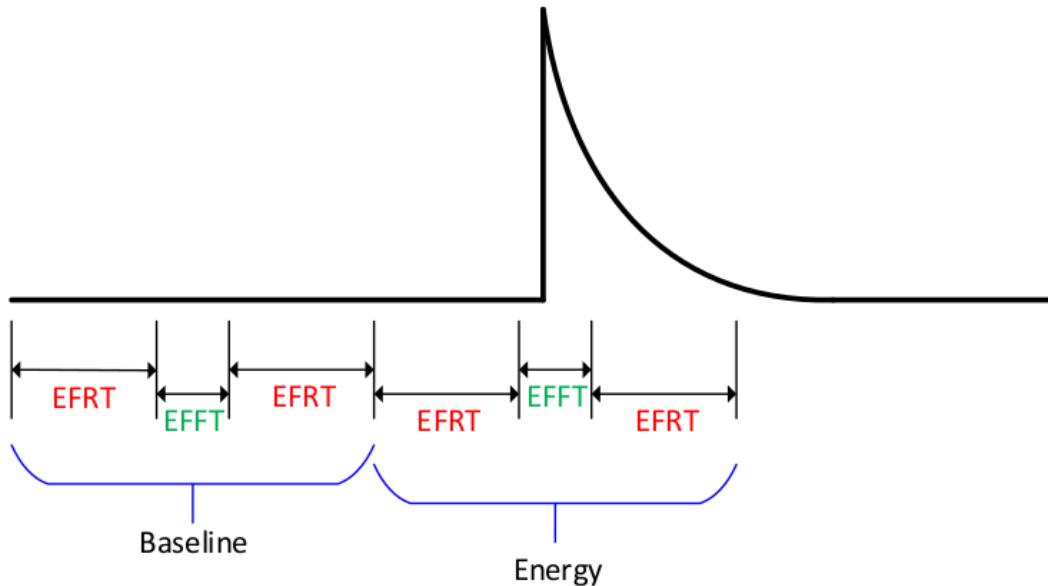
In the figure, EFRT stands for Energy Filter Rise Time and EFFT stands for Energy Filter Flat Top.

To compute energy filter response offline, the ideal settings are:

- Total trace length > $2 \times (2 \times \text{EFRT} + \text{EFFT})$
- Pre-trigger trace (Trace-delay) length > $(3 \times \text{EFRT} + \text{EFFT})$

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
the ADC trace display also includes the option to view a FFT of the acquired trace. This is useful to diagnose noise contributions.

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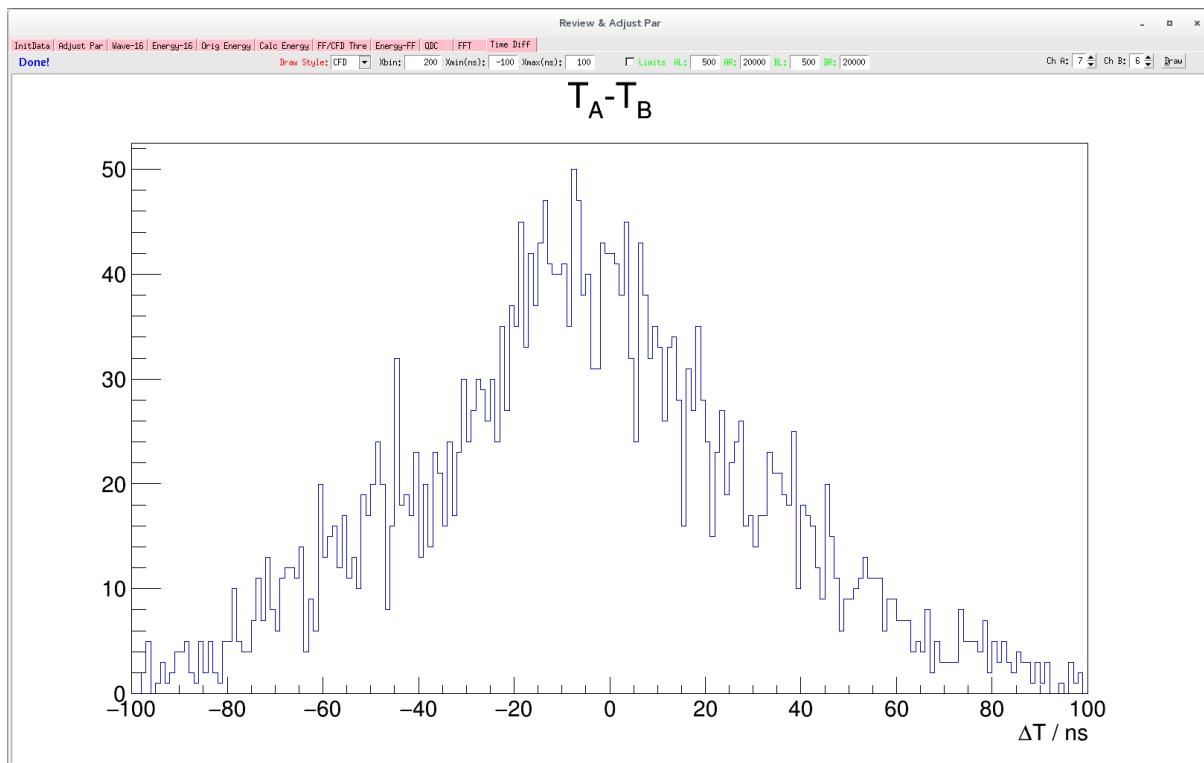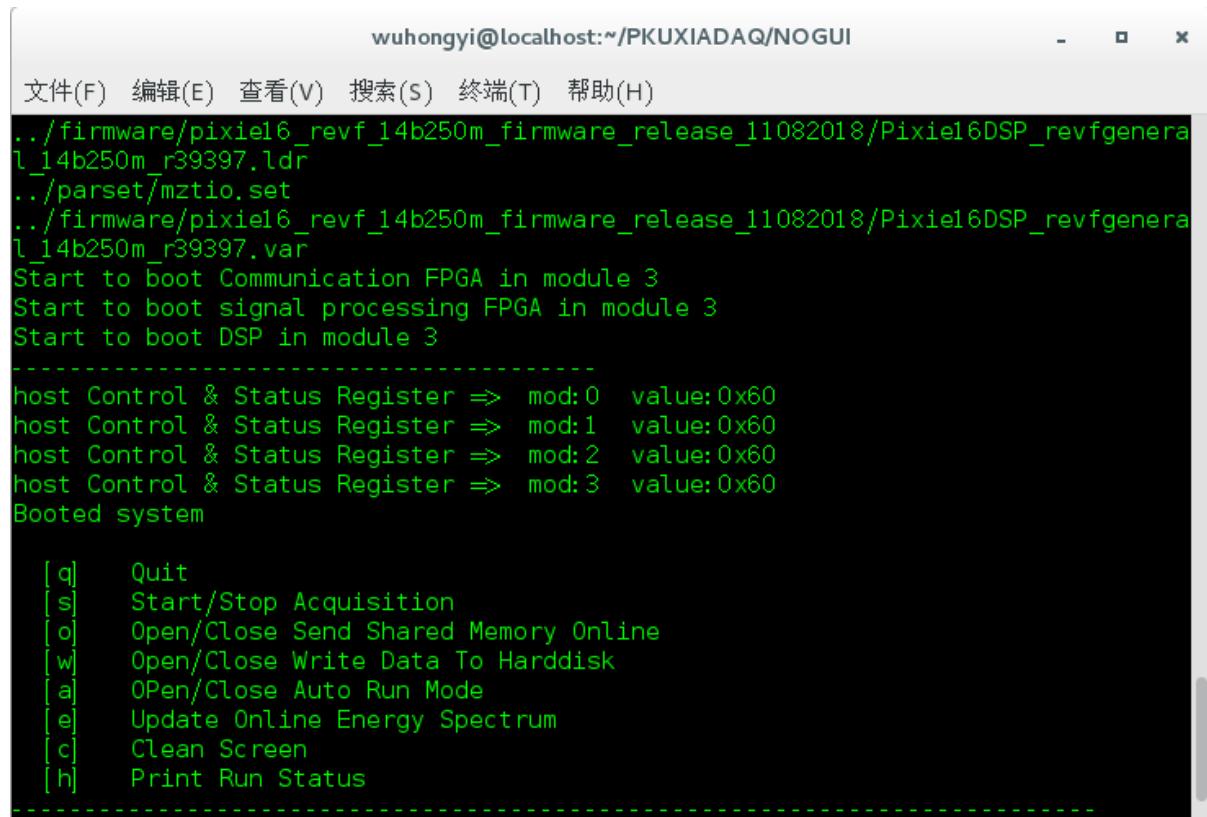
```
[q] Quit
[s] Start/Stop Acquisition
[o] Open/Close Send Shared Memory Online
[w] Open/Close Write Data To Harddisk
[a] Open/Close Auto Run Mode
[e] Update Online Energy Spectrum
[c] Clean Screen
[h] Print Run Status
```

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```
s
created the directory /home/wuhongyi/data/0320.
open: /home/wuhongyi/data/0320/data_R0320_M00.bin
open: /home/wuhongyi/data/0320/data_R0320_M01.bin
open: /home/wuhongyi/data/0320/data_R0320_M02.bin
open: /home/wuhongyi/data/0320/data_R0320_M03.bin
RUN START
SHM Open!
Running No. 320
Auto Run Mode: CLOSE
Send Shared Memory Online: OPEN
Write Data To Harddisk: OPEN
a
The times for each run: 180 s
h

[ q]    Quit
[ s]    Start/Stop Acquisition
[ o]    Open/Close Send Shared Memory Online
[ w]    Open/Close Write Data To Harddisk
[ a]    Open/Close Auto Run Mode
[ e]    Update Online Energy Spectrum
[ c]    Clean Screen
[ h]    Print Run Status
-----
Running No. 320
Auto Run Mode: OPEN - 180 s per run
Send Shared Memory Online: OPEN
Write Data To Harddisk: OPEN
```

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```
# Only use in NOGUI, unit: second
AutoRunModeTimes    180
```

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GDDAQ, 发布 beta

GDDAQ Online																					
File		Setup																			
CountRate			Alert			EnergyMonitor			File Path:			File Name:			Complete			RunStop		n0520 M11	
File Path:	/home/wuhongyi/data/	File Name:	data	Complete	RunStop	n0520	M11														
Monitor																					
ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	
00	0	0	00	0	0	00	0	0	00	0	0	00	0	0	00	0	0	00	0	0	00
01	0	0	01	0	0	01	0	0	01	0	0	01	0	0	01	0	0	01	0	0	01
02	0	0	02	0	0	02	0	0	02	0	0	02	0	0	02	0	0	02	0	0	02
03	0	0	03	0	0	03	0	0	03	0	0	03	0	0	03	0	0	03	0	0	03
04	0	0	04	0	0	04	0	0	04	0	0	04	0	0	04	0	0	04	0	0	04
05	0	0	05	0	0	05	0	0	05	0	0	05	0	0	05	0	0	05	0	0	05
06	0	0	06	0	0	06	0	0	06	0	0	06	0	0	06	0	0	06	0	0	06
07	0	0	07	0	0	07	0	0	07	0	0	07	0	0	07	0	0	07	0	0	07
08	0	0	08	0	0	08	0	0	08	0	0	08	0	0	08	0	0	08	0	0	08
09	0	0	09	0	0	09	0	0	09	0	0	09	0	0	09	0	0	09	0	0	09
10	0	0	10	0	0	10	0	0	10	0	0	10	0	0	10	0	0	10	0	0	10
11	0	0	11	0	0	11	0	0	11	0	0	11	0	0	11	0	0	11	0	0	11
12	0	0	12	0	0	12	0	0	12	0	0	12	0	0	12	0	0	12	0	0	12
13	0	0	13	0	0	13	0	0	13	0	0	13	0	0	13	0	0	13	0	0	13
14	0	0	14	0	0	14	0	0	14	0	0	14	0	0	14	0	0	14	0	0	14
15	0	0	15	0	0	15	0	0	15	0	0	15	0	0	15	0	0	15	0	0	15
ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	ch #	InRate/s	OutRate/s	File	Size/MB	ADC/MHz	
00	0	0	00	0	0	00	0	0	00	0	0	00	0	0	00	0	0	M00	0	100	
01	0	0	01	0	0	01	0	0	01	0	0	01	0	0	01	0	0	M01	0	100	
02	0	0	02	0	0	02	0	0	02	0	0	02	0	0	02	0	0	M02	0	100	
03	0	0	03	0	0	03	0	0	03	0	0	03	0	0	03	0	0	M03	0	100	
04	0	0	04	0	0	04	0	0	04	0	0	04	0	0	04	0	0	M04	0	100	
05	0	0	05	0	0	05	0	0	05	0	0	05	0	0	05	0	0	M05	0	100	
06	0	0	06	0	0	06	0	0	06	0	0	06	0	0	06	0	0	M06	0	100	
07	0	0	07	0	0	07	0	0	07	0	0	07	0	0	07	0	0	M07	0	100	
08	0	0	08	0	0	08	0	0	08	0	0	08	0	0	08	0	0	M08	0	100	
09	0	0	09	0	0	09	0	0	09	0	0	09	0	0	09	0	0	M09	0	100	
10	0	0	10	0	0	10	0	0	10	0	0	10	0	0	10	0	0	M10	0	100	
11	0	0	11	0	0	11	0	0	11	0	0	11	0	0	11	0	0	M11			
12	0	0	12	0	0	12	0	0	12	0	0	12	0	0	12	0	0	M12			
13	0	0	13	0	0	13	0	0	13	0	0	13	0	0	13	0	0				
14	0	0	14	0	0	14	0	0	14	0	0	14	0	0	14	0	0				
15	0	0	15	0	0	15	0	0	15	0	0	15	0	0	15	0	0				

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CHAPTER 11

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On the front panel of each Pixie-16 module, there are 16 analog signal input connectors, one LVDS I/O port, five digital I/O connectors as well as three LEDs near the bottom of the front panel. In addition, a sticker showing Pixie-16 model number (e.g., P16L-250-14, meaning the 14-bit, 250 MHz variant of the Pixie-16) is affixed to the top handle of the front panel, and another sticker indicating the serial number of the Pixie-16 module (e.g., S/N 1100) is placed at the bottom handle of the front panel.


11.1 Analog Signal Input Connectors(all revisions)

- **Connector Labels**
 - 0 to 15 for 16 channels
- **Connector Type**
 - SMB Jack

Each Pixie-16 module accepts 16 analog input signals, and each input connector is a SMB Jack (male contact) connector.

11.2 LVDS I/O Port (all revisions)






Connected Signals	4 LVDS pairs (Fo_{1p}/Fo_{1n} , Fi_{1p}/Fi_{1n} , Fi_{5p}/Fi_{5n} , Fo_{5p}/Fo_{5n} , see below for pin layout)
Signal Direction	Input or Output, software configurable
Cable Type	Cat 5 or Cat 6 (the same ones used for Ethernet)

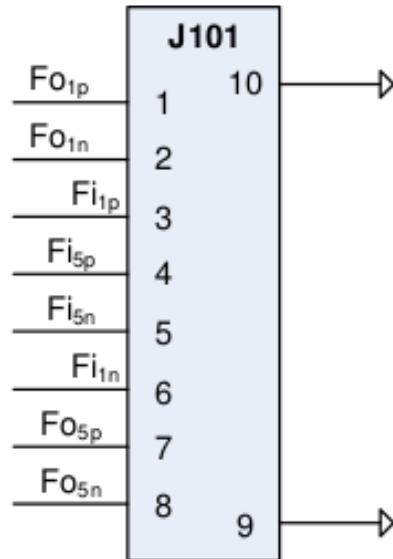Each Pixie-16 module is equipped with one LVDS I/O port on its front panel. LVDS stands for low voltage differential signaling. The LVDS I/O connector is a RJ45 connector, which implies that the same Cat 5 or Cat 6 Ethernet cables can be used to connect signals to or from this I/O port. However, no Ethernet connectivity is available through this Pixie-16 I/O port.

Four differential signal pairs, i.e., between pin-pairs Fo_{1p} / Fo_{1n} , Fi_{1p} / Fi_{1n} , Fi_{5p} / Fi_{5n} , and Fo_{5p} / Fo_{5n} , are available from this I/O port. Each pair can be configured as either an input or output signal.

11.3 Digital I/O Connectors(Rev. F only)

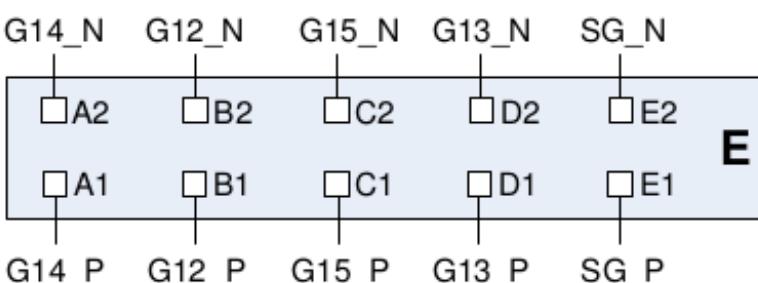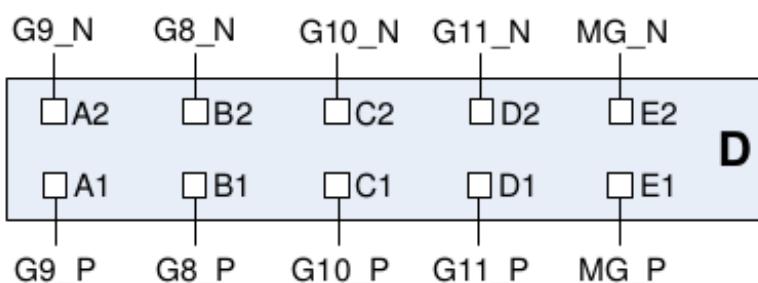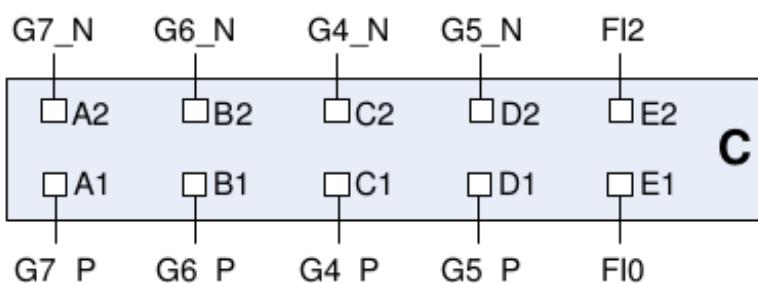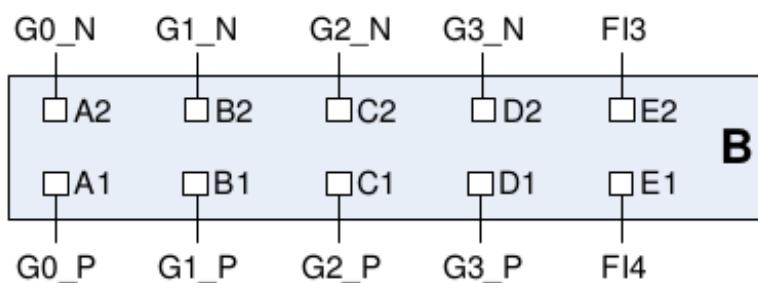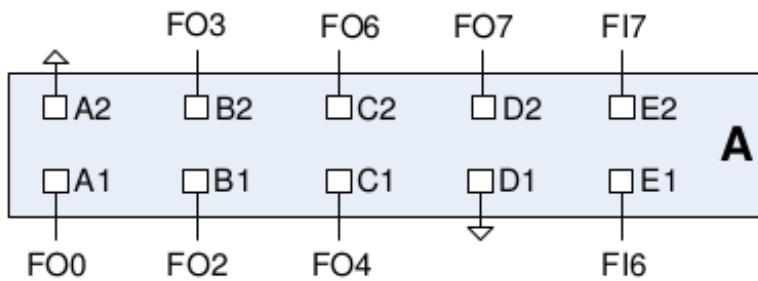The Pixie-16 Rev. F modules are equipped with five har-link[®] connectors on its front panel which act as digital I/O connectors. The 2mm pitch har-link[®] connector from HARTING is designed for high speed data transfer with

rates up to 2 Gbit/s. Its EMI shielding, shown in Figure 2mm pitch har-link® connector from HARTING, guarantees excellent performance in EM-polluted environment.

Each har-link® connector has 2 rows with 5 pins on each row, and is labelled using one of the five letters in red color font, from A to E. The signals connected to each pin of these five connectors are shown in Table *Rev. F Module's Digital I/O Connectors*.

Connector Type	har-link® (HARTING, 2mm pin spacing)
FI ₀ , FI ₂ , FI ₃ , FI ₄ , FI ₆ , FI ₇	TTL digital input signals (max. 5V)
FO ₀ , FO ₂ , FO ₃ , FO ₄ , FO ₆ , FO ₇	Digital outputs for test/debug purpose (TTL 5V)
Gx_P/Gx_N (x=0-15)	Channel Gate Inputs (0-15 for 16 channels) (LVDS format)
MG_P/MG_N	Module Gate Input (LVDS format)
SG_P/SG_N	Spare Gate Input (LVDS format)

Among them, FI₀, FI₂, FI₃, FI₄, FI₆, FI₇ are six TTL digital input signals. They can be signals like global fast trigger, global validation trigger, external clock, run inhibit, etc. The specific usage of each input pin is determined by the specific firmware that is downloaded to the Pixie-16 module (see Table 1-9 for input signals supported by the standard firmware). The six digital output signals, FO₀, FO₂, FO₃, FO₄, FO₆, FO₇, which are connected to six test output pins on the System FPGA of the Pixie-16, can be used to assist a user in the process of system setup. These test pins are connected to various internal signals of the Pixie-16 to provide insight of the current status of the system.

The Channel Gate Inputs (0-15 for 16 channels) are LVDS format input signals which independently gate the data acquisition of each of the 16 channels of a Pixie-16 module.

The Channel Gate signal is level sensitive signal, i.e., when the level of the Channel Gate Signal is logic high(1), the gate signal is effective; when the level of the Channel Gate Signal is logic low(0), the gate signal is not in use. In normal cases, the Channel Gate Signal is set up to veto the data acquisition in a given channel, i.e., at the time of the arrival of fast trigger in that channel, if the Channel Gate Signal is logic high(1), that fast trigger is discarded since it is vetoed. However, this type of logic can be reversed through setting corresponding registers in the FPGA via software. In such cases, the Channel Gate Signal is set up to validate the data acquisition in a given channel, i.e.,


at the time of fast trigger in that channel, only if the Channel Gate Signal is logic high(1) will that fast trigger be accepted to have the event recorded.

The Module Gate Input is a LVDS format signal that gates the data acquisition in all 16 channels of a Pixie-16 module. It is also a level sensitive signal, i.e., when the level of the Module Gate Signal is logic high(1), the gate signal is effective; when the level of the Module Gate Signal is logic low(0), the gate signal is not in use. In normal cases, the Module Gate Signal is set up to veto the data acquisition in all 16 channels, i.e., at the time of the arrival of fast trigger in any of the 16 channels, if the Module Gate Signal is logic high(1), that fast trigger of that channel is discarded since it is vetoed. However, this type of logic can be reversed through setting corresponding registers in the FPGA via software.

In such cases, the Module Gate Signal is set up to validate the data acquisition in all 16 channels, i.e., at the time of fast trigger in any of the 16 channel, only if the Module Gate Signal is logic high(1) will that fast trigger of that channel be accepted to have the event recorded.

The Spare Gate Input is a LVDS format signal that is reserved for special applications.

Such applications typically require development of custom firmware to support special functionalities of the Pixie-16 system.

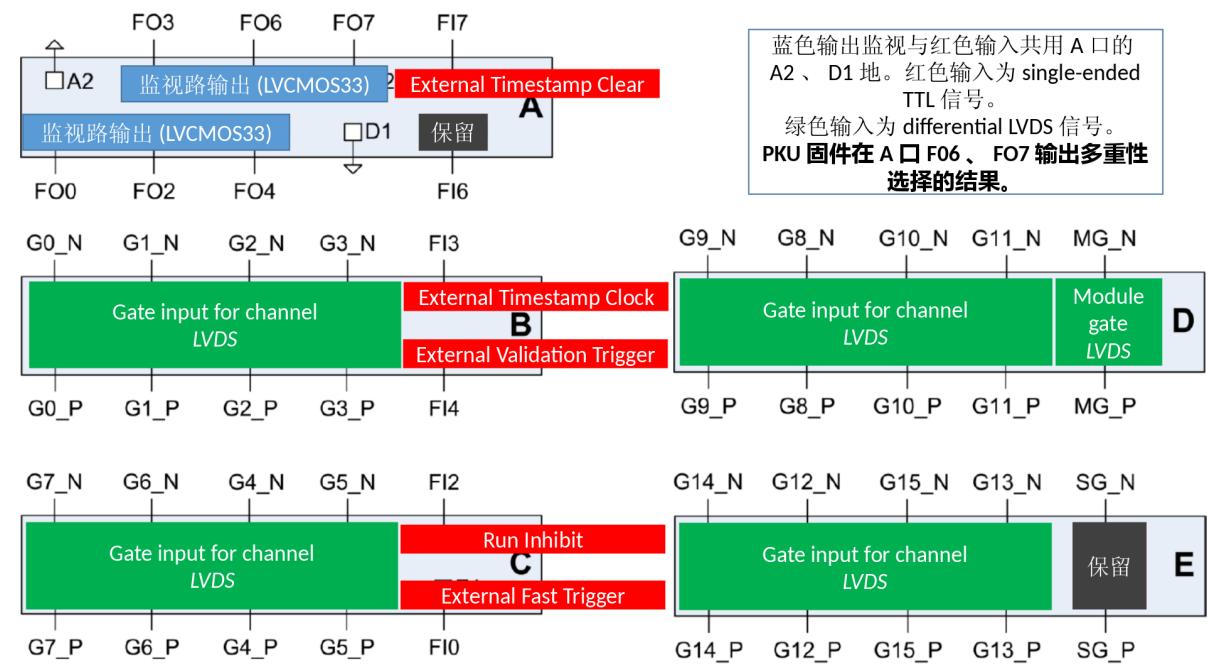11.4 Front Panel LEDs (all revisions)

Near the bottom of the Pixie-16 front panel, there are three LEDs. They are labelled as RUN, I/O, and ERR, respectively, from left to right. They correspond to three different colors, **green, yellow, and red**, respectively.

LED Name	Color	Function
RUN	Green	ON when run is in progress, and OFF if run is stopped or not started yet
I/O	Yellow	Flashing when there is I/O activity on the PCI bus between the Pixie-16 module and host computer
ERR	Red	ON when there is no more space in the External FIFO for storage of list mode event data, and OFF when there is sufficient space to store at least one more list mode event data (ON does not indicate any actual error condition. Rather, it simply indicates the External FIFO's FULL condition)


The **RUN LED** will be turned on when a run in the Pixie-16 module is in progress, and will be turned off when the run is stopped or not started yet.

The **I/O LED** will blink when there is I/O activity on the PCI bus between the Pixie-16 module and host computer.

The **ERR LED** is, in fact, not to signal any error condition in the Pixie-16 module. Instead, it is used to indicate whether or not the External FIFO of the Pixie-16 module is full. It will be ON when there is no more space in the External FIFO for storage of list mode event data, and OFF when there is sufficient space to store at least one more list mode event data. When the External FIFO is full, no more list mode event data can be written into it until the host software reads out part of the data in the External FIFO through the PCI bus.

11.5 3.3V I/O Connector(Rev. D only)

Connector Type	Single-ended, 2mm pin spacing
FO ₀ , FO ₂ , FO ₃ , FO ₄ , FO ₆ , FO ₇	Digital outputs for test/debug purpose (3.3V)

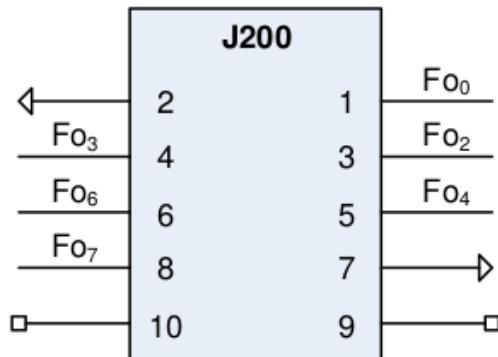
On Rev. D Pixie-16 modules, between analog input SMB connectors for channel 7 and channel 8, respectively, is the 3.3V I/O Connector(J200). It has 10 single-ended pins with 2mm spacing. Pins #1, 3, 4, 5, 6, and 8 are connected to six digital output signals from the System FPGA of the Pixie-16 module, i.e. FO0, FO2, FO3, FO4, FO6, FO7, mainly for the purpose of testing and debugging. Pins #2 and 7 are ground pins, and pins #9 and 10 are not in use.

11.6 GATE Inputs(Rev. D only)

Connector Type	Amphenol FCI® 55 Position Header, 2mm pin spacing
FI ₀ , FI ₂ , FI ₃ , FI ₄ , FI ₆ , FI ₇	TTL digital input signals (max. 5V)
Gx _{in+} /Gx _{in-} (x=0-15)	Channel Gate Inputs (0-15 for 16 channels) (LVDS format)
MG _{in+} /MG _{in-}	Module Gate Input (LVDS format)
SG _{in+} /SG _{in-}	Spare Gate Input (LVDS format)



On Rev. D Pixie-16 modules, between analog input SMB connectors for channel 11 and channel 12, respectively, is the GATE INPUTS connector. This connector is an Amphenol FCI 55 Position Header with 2mm pin spacing. The layout of these 55 pins is shown in Figure *Rev. D Module's GATE INPUTS Connector*. The 11 pins from the middle pin column (J150C) are all tied to the Ground. Among the first 8 rows of the GATE INPUTS connector, each differential pair of pins from the A/B columns (J150A/J150B) or the D/E columns (J150D/J150E) corresponds to one channel's GATE INPUT, which has the LVDS format, e.g. Gx_{in+}/Gx_{in-}(x=0-15). Differential pair of pins at J150A3/J150B3 is the Module Gate Input signal, MG_{in+}/MG_{in-}. Channel Gate Input signal can be used to veto or validate that given channel's own trigger. Module Gate Input signal works on the whole module level, i.e. it can be used to veto or validate all 16 channels' own trigger of that given module. Differential pair of pins at J150D3/J150E3 is the Spare Gate Input signal, SG_{in+}/SG_{in-}. Spare Gate Input signal can be used for special applications which require a custom firmware.

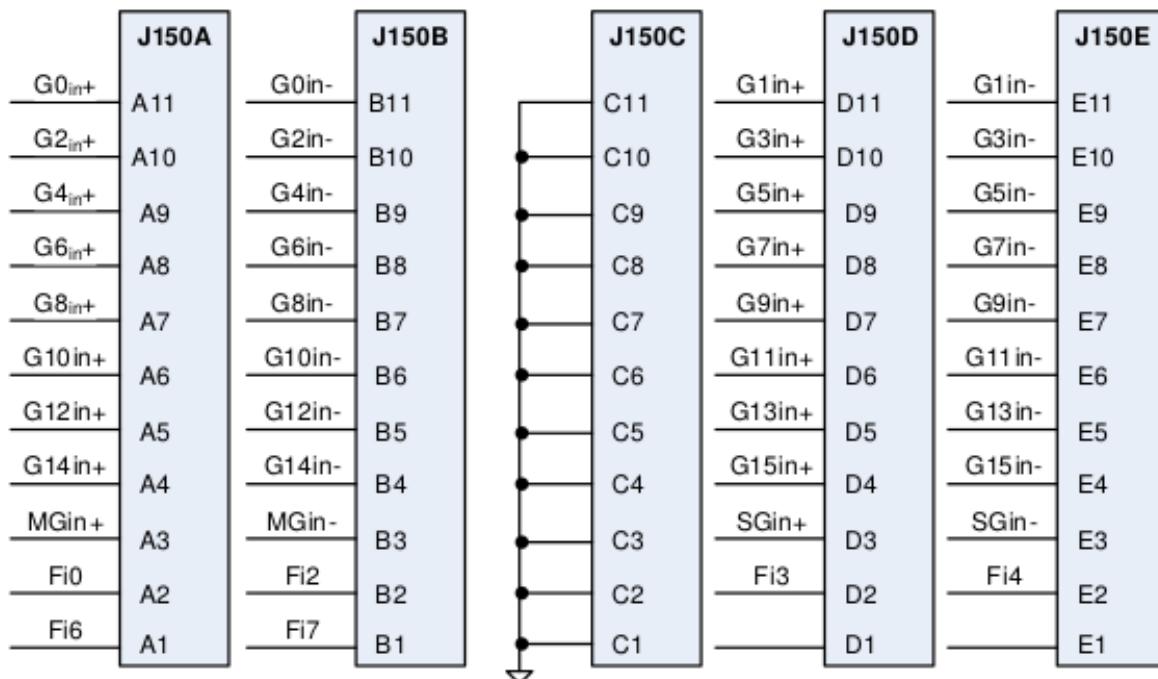On Rev. D Pixie-16 modules, the TTL digital input signals (max. 5V), i.e. FI₀, FI₂, FI₃, FI₄, FI₆, FI₇, are distributed among the bottom two rows of the GATE INPUTS Connector, as illustrated in Figure *Rev. D Module's GATE INPUTS Connector*.

11.7 3.3V I/O Connector(Rev. B and C only)

Connector Type	Single-ended, 2mm pin spacing
FI ₀ , FI ₂ , FI ₃ , FI ₄ , FI ₆ , FI ₇	TTL digital input signals (max. 5V)
FO ₀ , FO ₂ , FO ₃ , FO ₄ , FO ₆ , FO ₇	Digital outputs for test/debug purpose (3.3V)

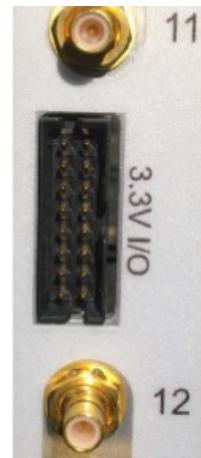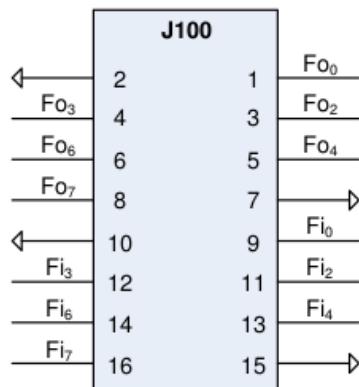On Rev. B and C Pixie-16 modules, between analog input SMB connectors for channel 11 and channel 12, respectively, is the 3.3V I/O Connector(J100). It has 16 single-ended pins with 2mm spacing. Pins #1, 3, 4, 5, 6, and 8 are connected to six digital output signals from the System FPGA of the Pixie-16 module, i.e. FO0, FO2, FO3, FO4, FO6, FO7, mainly for the purpose of testing and debugging. Pins #2, 7, 10 and 15 are ground pins. Pins #9, 11, 12, 13, 14 and 16 are connected to the six TTL digital input signals (max. 5V), i.e. FI0, FI2, FI3, FI4, FI6, FI7.

11.8 Digital Signals in Standard Firmware(all revisions)

The standard firmware of the Pixie-16 supports input and output of digital signals through its front panel I/O connectors, which were discussed earlier.

TTL digital input signals	Connected signals in standard firmware	Direction	Description
FI ₀	EXT_FASTTRIG	Input	External fast trigger signal
FI ₂	INHIBIT	Input	Run inhibit signal
FI ₃	EXT_TS_CLK	Input	External timestamp clock signal
FI ₄	EXT_VALIDTRIG	Input	External validation signal
FI ₆	not used		
FI ₇	EXT_TS_CLR	Input	External timestamp clear signal

Table *TTL Digital Input Signals* shows the five TTL digital input signals supported by the Pixie-16 standard firmware.

Among them, the signals **EXT_TS_CLK** and **EXT_TS_CLR** are used for external timestamping in the Pixie-16, i.e. the Pixie-16 accepting an external clock signal(the frequency of this external clock is not recommended to exceed about 20 MHz in order to avoid the clock signal integrity issue), counting such clock signal with a 48-bit counter, and outputting such counter value to the list mode data stream when an event trigger occurs.

The external timestamping is useful for synchronizing the Pixie-16 data acquisition system with another data acquisition system through correlating the external timestamps of the events recorded by both systems.

The **INHIBIT** signal is used by an external system to inhibit the data acquisition run in a Pixie-16 system when synchronization requirement is enabled in the Pixie-16 modules. It is a level sensitive signal, i.e. when the **INHIBIT**

signal is at the logic high level, the run in the Pixie-16 won't start. Only when the **INHIBIT** signal goes to the logic low level will the run start in the Pixie-16. During the run, if the **INHIBIT** signal returns to the logic high level, the run will be aborted.

The **EXT_FASTTRIG** signal is the external fast trigger signal, which can be used to replace the local fast trigger for recording events in the Pixie-16 modules. The **EXT_VALIDTRIG** signal is the external validation signal, which can be used to validate events in the Pixie-16 modules.

Connector J101 Pins	Connected signals in standard firmware	Direction	Description
F _{01p} /F _{01n}	not used		
F _{11p} /F _{11n}	LVDS_VALIDTRIG	Input	External validation trigger signal in LVDS format
F _{15p} /F _{15n}	LVDS_FASTTRIG	Input	External fast trigger signal in LVDS format
F _{05p} /F _{05n}	SYNC_LVDS_FP	Output	Pixie-16 synchronization output signal in LVDS format (to synchronize with other DAQ systems)

Table *Connector J101 LVDS I/O Port Signals* shows the Pixie-16 connector J101 LVDS I/O port signals. This J101 LVDS I/O port can use the regular Ethernet cable for connection but it does not have Ethernet connectivity. Among the four LVDS pairs available from this J101 port, one pair is currently not in use, two pairs are used for input and one pair is used for output. The **LVDS_VALIDTRIG** is the external validation trigger signal in LVDS format, and the **LVDS_FASTTRIG** is the external fast trigger signal in LVDS format. The **SYNC_LVDS_FP** is an output signal from the Pixie-16 module to indicate to external data acquisition systems the synchronization status of the Pixie-16 system so that both data acquisition systems can be synchronized.

TTL digital output signals	Connected signals in standard firmware		Direction	Description	
F _{O0}	FTRIG_DELAY	FTRIG_DELAY	Output	Delayed local fast trigger of one of the 16 channels	Delayed local fast trigger of one of the 16 channels
F _{O2}	FTRIG_VAL	VETO_CE	Output	Validated, delayed local fast trigger one of the 16 channels	Stretched veto trigger of one of the 16 channels
F _{O3}	ETRIG_CE	LDPMFULL	Output	Stretched external global validation trigger of one of the 16 channels	Module level dual port memory (DPM) full status flag
F _{O4}	CHANTRIG_CE	SDPMFULL	Output	Stretched channel validation trigger of one of the 16 channels	System level dual port memory (DPM) full status flag
F _{O6}	FTIN_OR	FTIN_OR	Output	OR of 16 local fast triggers	OR of 16 local fast triggers
F _{O7}	TEST_SEL	TEST_SEL	Output	Selected test signal	Selected test signal

Table *TTL Digital Output Signals* lists the six Pixie-16 TTL digital output signals. Two groups of six output signals can be chosen through software settings (see bits [14:12] and [19:16] of TrigConfig0). The last output signal TEST_SEL

can be further selected through software settings. More details about these signals will be provided in later sections of this manual.

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


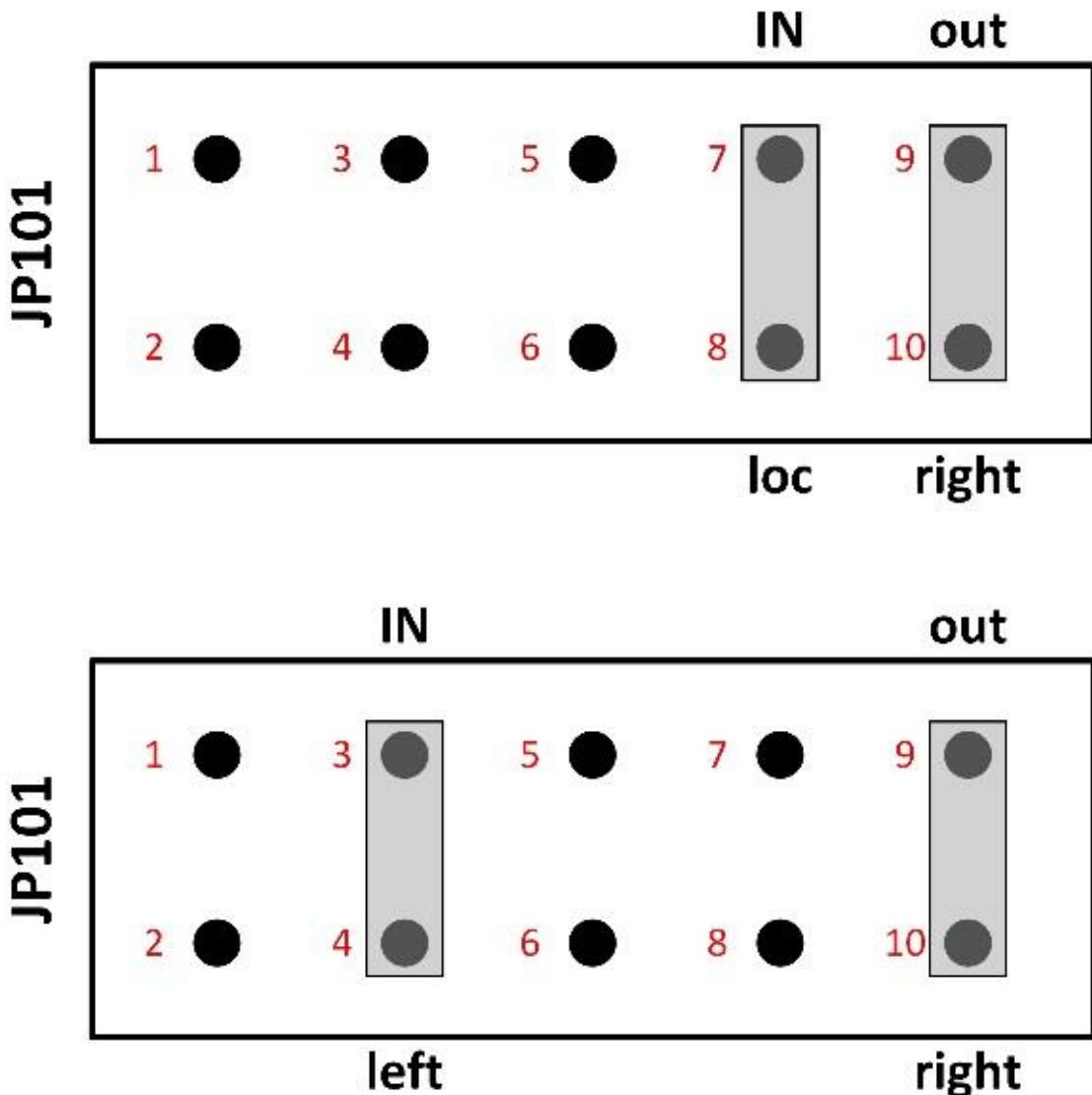
Top: Daisy-chain clock master (leftmost module)
 Bottom: Daisy-chain clock repeater

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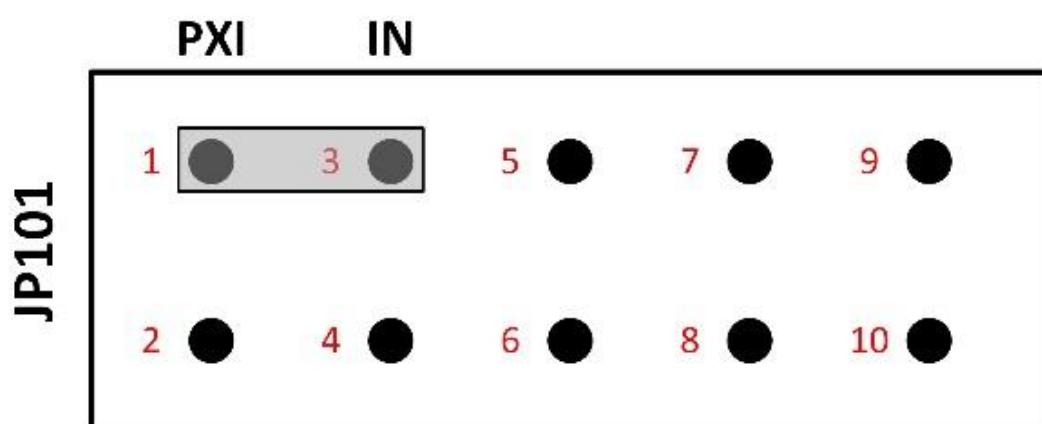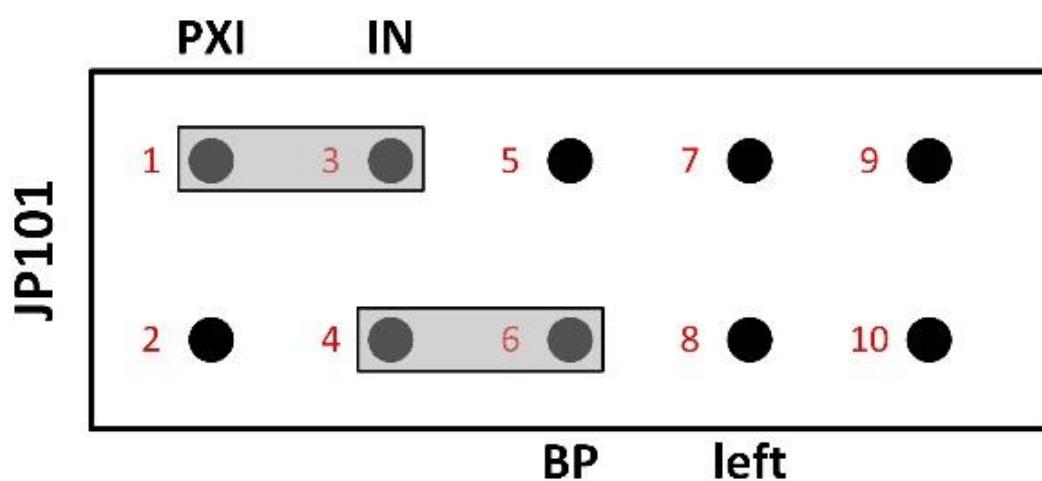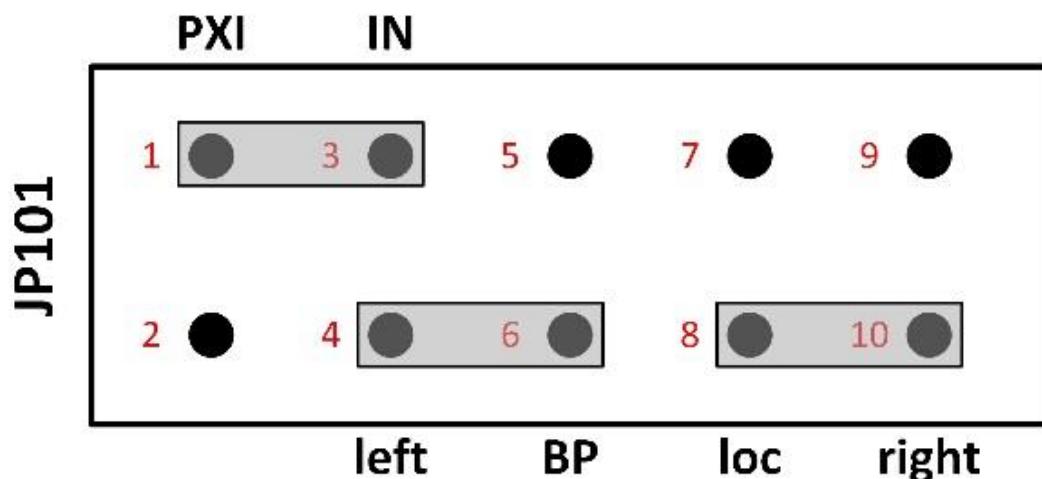System Director Module	Connect pins 1 and 3, 4 and 6, 8 and 10.
Crate Master Module	Connect pins 1 and 3, 4 and 6.
Regular Module	Connect pins 1 and 3.

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JP1	Connect pins 1 and 2 for “P16”
JP20	Connect pins 2 and 3, 6 and 7
JP40	Connect pins 2 and 3, 6 and 7
JP60	Connect pins 1 and 2, 7 and 8
JP21	Connect pins 2 and 3
JP41	Connect pins 2 and 3
JP61	Connect pins 1 and 2
JP100	Connect pins 2 and 3 (connect to J4)
JP101	Connect pins 2 and 3 (connect to J4)
JP102	Connect pins 2 and 3 (connect to J4)
JP103	Connect pins 2 and 3 (connect to J4)
JP104	Connect pins 2 and 3 (connect to J4)
JP105	Connect pins 2 and 3 (connect to J4)

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JP1	Connect pins 2 and 3 for “loc”
JP20	Connect pins 1 and 2, 7 and 8
JP40	Connect pins 1 and 2, 7 and 8
JP60	Connect pins 2 and 3, 6 and 7
JP21	Don’t connect any pin
JP41	Don’t connect any pin
JP61	Connect pins 1 and 2
JP100	Connect pins 2 and 3 (connect to J4)
JP101	Connect pins 2 and 3 (connect to J4)
JP102	Connect pins 2 and 3 (connect to J4)
JP103	Connect pins 2 and 3 (connect to J4)
JP104	Connect pins 2 and 3 (connect to J4)
JP105	Connect pins 2 and 3 (connect to J4)

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
14.1 6 pin JTAG connector

You should now try to locate the 6-pin vertical Pixie-16 CPLD JTAG header J8 (circled in red in the screenshot below).

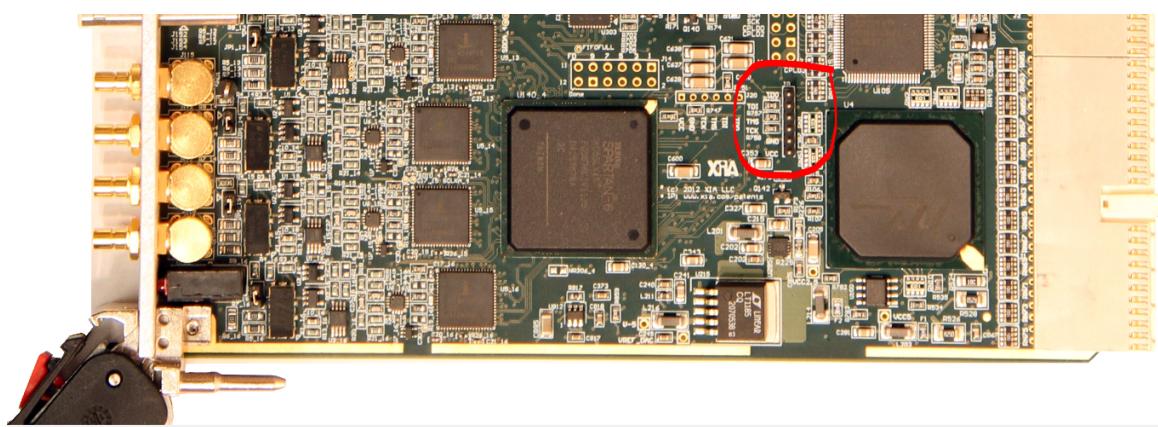The next step is to identify each of the 6 pins of the Pixie-16 CPLD JTAG header J8. When facing the top side of the Pixie-16 board, the 6 pins of the J8 connector can be identified as follows from top to bottom: **TDO**, **TDI**, **TMS**, **TCK**, **GND**, and **VCC**. It is very important to match these pins with the colored cables from the green PCB connector that is plugged into the Xilinx Platform USB programmer.

One issue here is that the Pixie-16 CPLD JTAG header J8 has 2 mm pin spacing, and the colored cables from the green PCB connector of the Xilinx programmer might have too large header to be plugged into the J8 pins directly. In that case, an adapter might be needed to convert those colored cable headers to 2 mm pin spacing first.

14.2 update the CPLD

Turn on your Pixie16 crate and the LED on the Xilinx Platform Cable USB should turn green. If the LED does not turn green check your cable connection. The grey INIT cable is purposely left unconnected. Your start menu should now include the folder Xilinx ISE Design Suite 14.7. Run the 32 bit or 64 bit version of “iMPACT”. Click on **File->Initialize Chain** and the program should find the **xc2c256 CPLD**. Click on “**No**” in the window “**Auto Assign Configuration Files Query**”, then click on “**Cancel**”. Left click on the Xilinx part in the window once, then right click and select “**Assign New Configuration Files**” and **select the new “jed” file**. Left click on the Xilinx part, right click and then click on “**Program**”. A device programming properties window appears. Ensure that only “**Verify**” and “**Erase before Programming**” are checked, then click on “**OK**”. The program should now say “**program succeeded**” in blue at the bottom of the window. The CPLD is now reprogrammed and the Pixie16 is ready to be tested.

I would suggest testing out your first board before reprogramming the CPLDs on all of your boards. Once you are in the batch mode of reprogramming and you have the cable connected (the USB programmer LED is green), you only need to click on “**Program**” to reprogram CPLD on your next Pixie16.

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CHAPTER 22

XIA API

It from **Programmer's Manual Digital Gamma Finder (DGF) PIXIE-16 Version 1.40, October 2009**

```
// Configure modules for communication in PXI chassis
// Use this function to configure the Pixie-16 modules in the PXI chassis.
// NumModules is the total number of Pixie-16 modules installed in the system.
// →PXISlotMap is the pointer to an array that must have at least as many entries as
// →there are Pixie-16 modules in the chassis.
// PXISlotMap serves as a simple mapping of the logical module number and the
// →physical slot number that the modules reside in. The logical module number runs
// →from 0. For instance, in a system with 5 Pixie-16 modules, these 5 modules may
// →occupy slots 3 through 7. The user must fill PXISlotMap as follows: PXISlotMap =
// →{3, 4, 5, 6, 7 ...} since module number 0 resides in slot number 3, etc. To find
// →out in which slot a module is located, any piece of subsequent code can use the
// →expression PXISlotMap[ModNum], where ModNum is the logic module number.
// OfflineMode is used to indicate to the API whether the system is running in
// →OFFLINE mode (1) or ONLINE mode (0). OFFLINE mode is useful for situations where
// →no Pixie-16 modules are present but users can still test their calls to the API
// →functions in their application software.
// This function must be called as the first step in the boot process. It makes
// →the modules known to the system and "opens" each module for communication.
// The function relies on an initialization file (pxisys.ini) that contains
// →information about the Host PC's PCI buses, including the slot enumeration scheme.
// →XIA's software distribution normally puts this file under the same folder as
// →Pixie-16 software installation folder. However, the user has the flexibility of
// →putting it in other folders by simply changing the definition of the string
// →PCISysIniFile_Windows or PCISysIniFile_Linux in the header part of the file
// →pixie16sys.c, depending on which operating system is being used.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16InitSystem (
    unsigned short NumModules,      // total number of Pixie16 modules in the
    //system
    unsigned short *PXISlotMap,     // an array containing the PXI slot number for
    //each pixie16 module
    unsigned short OfflineMode ); // specify if the system is in offline mode
```

```
// Release user virtual addressees assigned to modules
// Use this function to release the user virtual addressees that are assigned to
// →Pixie-16 modules when these modules are initialized by function
// →Pixie16InitSystem. This function should be called before a user's application
// →exits.
```

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```
// If ModNum is set to less than the total number of modules in the system, only
// the module specified by ModNum will be closed. But if ModNum is equal to the
// total number of modules in the system, e.g. there are 5 modules in the chassis
// and ModNum = 5, then all modules in the system will be closed altogether. Note
// that the modules are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ExitSystem (
    unsigned short ModNum ); // module number
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadModuleInfo (
    unsigned short ModNum, // module number
    unsigned short *ModRev, // returned module revision
    unsigned int *ModSerNum, // returned module serial number
    unsigned short *ModADCBits, // returned module ADC bits
    unsigned short *ModADCMSPS ); // returned module ADC sampling rate
```

```
// Boot modules so that they can be set up for data taking
// Use this function to boot Pixie-16 modules so that they can be set up for data
// taking. The function downloads to the Pixie-16 modules the communication FPGA
// configurations, signal processing FPGA configurations, trigger FPGA
// configurations (Revision A modules only), executable code for the digital signal
// processor (DSP), and DSP parameters.
// The FPGA configurations consist of a fixed number of words depending on the
// hardware mounted on the modules; the DSP codes have a length which depends on
// the actual compiled code; and the set of DSP parameters always consists of 1280
// 32-bit words for each module. The host software has to make the names of those
// boot data files on the hard disk available to the boot function.
// If ModNum is set to be less than the total number of modules in the system,
// only the module specified by ModNum will be booted. But if ModNum is equal to
// the total number of modules in the system, e.g. there are 5 modules in the
// chassis and ModNum = 5, then all modules in the system will be booted.
// The boot pattern is a bit mask (shown below) indicating which on-board chip
// will be booted. Under normal circumstances, all on-board chips should be booted,
// i.e. the boot pattern would be 0x7F. For Rev-B, C, D modules, bit 1, i.e., "Boot
// trigger FPGA", will be ignored even if that bit is set to 1.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BootModule (
    char *ComFPGAConfigFile, // name of communications FPGA configuration
    file
    char *SPFPGAConfigFile, // name of signal processing FPGA
    configuration file
    char *TrigFPGAConfigFile, // name of trigger FPGA configuration file
    char *DSPCodeFile, // name of executable code file for digital
    signal processor (DSP)
    char *DSPParFile, // name of DSP parameter file
    char *DSPVarFile, // name of DSP variable names file
    unsigned short ModNum, // pixie module number
    unsigned short BootPattern ); // boot pattern bit mask
```

```
// Acquire ADC traces in single or multiple modules
// Use this function to acquire ADC traces from Pixie-16 modules. Specify the
// module using ModNum. If ModNum is set to be less than the total number of
// modules in the system, only the module specified by ModNum will have its ADC
// traces acquired. But if ModNum is equal to the total number of modules in the
// system, then all modules in the system will have their ADC traces acquired.
// After the successful return of this function, the DSP's internal memory will be
// filled with ADC trace data. A user's application software should then call
// another function Pixie16ReadSglChanADCTrace to read the ADC trace data out to
// the host computer, channel by channel.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireADCTrace (
    unsigned short ModNum ); // module number
```

```

// Read ADC trace data from a channel in a module
// Use this function to read ADC trace data from a Pixie-16 module. Before calling
// this function, another function Pixie16AcquireADCTrace should be called to fill
// the DSP internal memory first. Also, the host code should allocate appropriate
// amount of memory to store the trace data. The ADC trace data length for each
// channel is 8192. Since the trace data are 16-bit unsigned integers (actually
// only the lower 14-bit contains real data due to the on-board 14-bit ADC), two
// consecutive 16-bit words are packed into one 32-bit word in the DSP internal
// memory. So for each channel, 4096 32-bit words are read out first from the DSP,
// and then each 32-bit word is unpacked to form two 16-bit words.
// Specify the module using ModNum and the channel on the module using ChanNum.
// Note that both the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanADCTrace (
    unsigned short *Trace_Buffer, // trace data
    unsigned int Trace_Length,   // trace length
    unsigned short ModNum,      // module number
    unsigned short ChanNum );   // channel number

```

```

// Transfer data between host and DSP internal memory
// Use this function to directly transfer data between the host and the DSP
// internal memory of a Pixie-16 module. The DSP internal memory is split into two
// blocks with address range 0x40000 to 0x4FFFF for the first block and address
// range 0x50000 to 0x5FFFF for the second block. Within the first block, address
// range 0x40000 to 0x49FFF is reserved for program memory and shouldn't be
// accessed directly by the host. Address range 0x4A000 to 0x4A4FF is used by the
// DSP I/O parameters which are stored in the configuration files (.set files) in
// the host. Within this range, 0x4A000 to 0x4A33F can be both read and written,
// but 0x4A340 to 0x4A4FF can only be read but not written. The remaining address
// range (0x4A500 to 4FFFF) in the first block and the entire second block (0x50000
// to 0x5FFFF) should only be read but not written by the host. Use Direction = 1
// for read and Direction = 0 for write.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16IMbufferIO (
    unsigned int *Buffer,           // buffer data
    unsigned int NumWords,          // number of buffer words to read or write
    unsigned int Address,           // buffer address
    unsigned short Direction,       // I/O direction
    unsigned short ModNum );        // module number

```

```

// Transfer data between host and DSP external memory
// Use this function to directly read data from or write data to the on-board
// external memory of a Pixie-16 module. The valid memory address is from 0x0 to
// 0x7FFFF (32-bit wide). Use Direction = 1 for read and Direction = 0 for write.
// The external memory is used to store the histogram data accumulated for each of
// the 16 channels of a Pixie-16 module. Each channel has a fixed histogram length
// of 32768 words (32-bit wide), and the placement of the histogram data in the
// memory is in the same order of the channel number, i.e. channel 0 occupies
// memory address 0x0 to 0x7FFF, channel 1 occupies 0x8000 to 0xFFFF, and so on.
// NOTE: another function Pixie16ReadHistogramFromModule can also be used to read
// out the histograms except that it needs to be called channel by channel.
// In Rev-A modules, part of the external memory is also used to store the list
// mode data in ping-pong buffering mode. This function can be used to read list
// mode data from the buffers.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EMbufferIO (
    unsigned int *Buffer,           // buffer data
    unsigned int NumWords,          // number of buffer words to read or write
    unsigned int Address,           // buffer address
    unsigned short Direction,       // I/O direction
    unsigned short ModNum );        // module number

```

```

// Start a list mode data acquisition run
// Use this function to start a list mode data acquisition run in Pixie-16 modules.
// List mode run is used to collect data on an event-by-event basis, gathering
// energies, timestamps, pulse shape analysis values, and waveforms, for each event.
// Runs will continue until a preset number of events are reached or the user
// terminates the run by calling function Pixie16EndRun. Once the run is progress, 145
// if the run is set to terminate after a given number of events have been
// accumulated, another function, Pixie16CheckRunStatus, should be called to check
// if the run has finished. To start the data acquisition this function has to be
// called for every Pixie-16 module in the system. If all modules are to run

```

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```

// Use mode=NEW_RUN (=1) to erase histograms and statistics information before_
// launching the new run. Note that this will cause a start up delay of up to 1_
// millisecond. Use mode=RESUME_RUN (=0) to resume an earlier run. This mode has a_
// start up delay of only a few microseconds.
// For Rev-A modules, currently there are 4 list mode run types supported. They_
// are 0x100 (general purpose run), 0x101 (without waveforms), 0x102 (without_
// auxiliary data) and 0x103 (energy and timestamp only).
// For Rev-B, C, D modules, there are only one list mode run type supported, that_
// is, 0x100. However, different output data options can be chosen by enabling or_
// disabling different CHANCSRA bits.
// Histograms and statistics data are updated incrementally from run to run_
// provided RESUME_RUN mode is used.
// ModNum is the module number which starts counting at 0. If ModNum is set to be_
// less than the total number of modules in the system, only the module specified_
// by ModNum will have its list mode run started. But if ModNum is set to equal to_
// the total number of modules in the system, then all modules in the system will_
// have their runs started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartListModeRun (
    unsigned short ModNum,           // module number
    unsigned short RunType,         // run type
    unsigned short mode );          // run mode

```

```

// Start a MCA histogram mode data acquisition run
// Use this function to begin a data acquisition run that accumulates energy_
// histograms, one for each channel. It launches a data acquisition run in which_
// only energy information is preserved and histogrammed locally to each channel.
// Call this function for each Pixie-16 module in the system. The last module_
// addressed will allow the actual data acquisition to begin. Histogram run can be_
// self-terminating when the elapsed run time exceeds the preset run time, or the_
// user can prematurely terminate the run by calling Pixie16EndRun. On completion,_
// final histogram and statistics data will be available.
// Use mode=NEW_RUN (=1) to erase histograms and statistics information before_
// launching the new run. Use mode=RESUME_RUN (=0) to resume an earlier run.
// ModNum is the module number which starts counting at 0. If ModNum is set to be_
// less than the total number of modules in the system, only the module specified_
// by ModNum will have its histogram run started. But if ModNum is set to be equal_
// to the total number of modules in the system, then all modules in the system_
// will have their runs started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartHistogramRun (
    unsigned short ModNum,           // module number
    unsigned short mode );          // run mode

```

```

// Check status of a data acquisition run
// Use this function to check the run status of a Pixie-16 module while a list_
// mode data acquisition run is in progress. If the run is still in progress_
// continue polling.
// If the return code of this function indicates the run has finished, there might_
// still be some data in the external memory (Rev-A modules) or external FIFO (Rev-
// B, C, D modules) that need to be read out to the host. In addition, final run_
// statistics and histogram data are available for reading out too.
// In MCA histogram run mode, this function can also be called to check if the run_
// is still in progress even though it is normally self-terminating.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckRunStatus (
    unsigned short ModNum );        // Pixie module number

```

```

// Stop a data acquisition run
// Use this function to end a histogram run, or to force the end of a list mode_
// run. In a multi-module system, if all modules are running synchronously, only_
// one module needs to be addressed this way. It will immediately stop the run in_
// all other module in the system.

```

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```
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EndRun (
    unsigned short ModNum ); // Pixie module number
```

```
// Compute input count rate
// Use this function to calculate the input count rate on one channel of a Pixie-
// →16 module. This function does not communicate with Pixie-16 modules. Before
// →calling this function, another function, Pixie16ReadStatisticsFromModule, should
// →be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer
// →words (32-bit). The *Statistics array is filled with data from a Pixie-16 module
// →after calling function Pixie16ReadStatisticsFromModule. ModNum is the module
// →number which starts counting at 0. ChanNum is the channel number which starts
// →counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeInputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );
```

```
// Compute output count rate of a channel
// Use this function to calculate the output count rate on one channel of a Pixie-
// →16 module. This function does not communicate with Pixie-16 modules. Before
// →calling this function, another function, Pixie16ReadStatisticsFromModule, should
// →be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer
// →words (32-bit). The *Statistics array is filled with data from a Pixie-16 module
// →after calling function Pixie16ReadStatisticsFromModule. ModNum is the module
// →number which starts counting at 0. ChanNum is the channel number which starts
// →counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeOutputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );
```

```
// Compute live time that a channel accumulated in a run
// Use this function to calculate the live time that one channel of a Pixie-16
// →module has spent on data acquisition. This function does not communicate with
// →Pixie-16 modules. Before calling this function, another function,
// →Pixie16ReadStatisticsFromModule, should be called to read statistics data from
// →the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer
// →words (32-bit). The *Statistics array is filled with data from a Pixie-16 module
// →after calling function Pixie16ReadStatisticsFromModule. ModNum is the module
// →number which starts counting at 0. ChanNum is the channel number which starts
// →counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeLiveTime (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );
```

```
// Compute number of events processed by a channel
// Use this function to calculate the number of events that have been processed by
// →a Pixie-16 module during a data acquisition run. This function is only used by
// →Rev-A modules. This function does not communicate with Pixie-16 modules. Before
// →calling this function, another function, Pixie16ReadStatisticsFromModule, should
// →be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer
// →words (32-bit). The *Statistics array is filled with data from a Pixie-16 module
// →after calling function Pixie16ReadStatisticsFromModule. ModNum is the module
// →number which starts counting at 0. ChanNum is the channel number which starts
// →counting at 0.
```

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```
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeProcessedEvents (
    unsigned int *Statistics,
    unsigned short ModNum );
```

```
// Compute real time that a channel accumulated in a run
// Use this function to calculate the real time that a Pixie-16 module has spent
// on data acquisition. This function does not communicate with Pixie-16 modules.
// Before calling this function, another function, Pixie16ReadStatisticsFromModule,
// should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer
// words (32-bit). The *Statistics array is filled with data from a Pixie-16 module
// after calling function Pixie16ReadStatisticsFromModule. ModNum is the module
// number which starts counting at 0. ChanNum is the channel number which starts
// counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeRealTime (
    unsigned int *Statistics,
    unsigned short ModNum );
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16complexFFT (
    double *data,
    unsigned int length );

// Test one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_TstBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_SetBit (
    unsigned short bit,
    unsigned short value );

// Clear one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_ClrBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_SetBit (
    unsigned short bit,
    unsigned int value );

// Clear one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_ClrBit (
    unsigned short bit,
    unsigned int value );

// Test one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_TstBit (
    unsigned short bit,
    unsigned int value );
```

```
// Program on-board DACs
// Use this function to reprogram the on-board digital to analog converters (DAC)
// of the Pixie-16 modules. In this operation the DSP uses data from the DSP
// parameters that were previously downloaded.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SetDACS (
    unsigned short ModNum );
```

```
// Program on-board signal processing FPGAs
// Use this function to program the on-board signal processing FPGAs of the Pixie-
// →16 modules. After the host computer has written the DSP parameters to the DSP_
// →memory, the DSP needs to write some of these parameters to the FPGAs. This_
// →function makes the DSP perform that action.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ProgramFippi (
    unsigned short ModNum );
```

```
// Adjust DC-offsets in single or multiple modules
// Use this function to adjust the DC-offsets of Pixie-16 modules. Specify the_
// →module using ModNum. If ModNum is set to be less than the total number of_
// →modules in the system, only the module specified by ModNum will have its DC-
// →offsets adjusted. But if ModNum is set to be equal to the total number of_
// →modules in the system, then all modules in the system will have their DC-offsets_
// →adjusted.
// After the DC-offset levels have been adjusted, the baseline level of the_
// →digitized input signals will be determined by the DSP parameter BaselinePercent._
// →For instance, if BaselinePercent is set to 10(%), the baseline level of the_
// →input signals will be ~ 1638 on the 14-bit ADC scale (minimum: 0; maximum:_
// →16383).
// The main purpose of this function is to ensure the input signals fall within_
// →the voltage range of the ADCs to ensure all input signals can be digitized by_
// →the ADCs properly.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AdjustOffsets (
    unsigned short ModNum );
```

```
// Acquire baselines from a module
// Use this function to acquire baselines from Pixie-16 modules. Specify the_
// →module using ModNum. If ModNum is set to be less than the total number of_
// →modules in the system, only the module specified by ModNum will have its_
// →baselines acquired. But if ModNum is set to be equal to the total number of_
// →modules in the system, then all modules in the system will have their baselines_
// →acquired.
// After the successful return of this function, the DSP's internal memory will be_
// →filled with baselines data. Users should then call another function_
// →Pixie16ReadSglChanBaselines to read the baselines data out to the host computer,_
// →channel by channel.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireBaselines (
    unsigned short ModNum );           // module number
```

```
// Read baselines from a channel in a module
// Use this function to read baselines data from a Pixie-16 module. Before calling_
// →this function, another function Pixie16AcquireBaselines should be called to fill_
// →the DSP internal memory first. Also, the host code should allocate appropriate_
// →amount of memory to store the baseline data. The baselines data length for each_
// →channel is 3640. In the DSP internal memory, each baseline data is a 32-bit IEEE_
// →floating point number. After being read out to the host, this function will_
// →convert each baseline data to a decimal number. In addition to baseline values,_
// →timestamps corresponding to each baseline were also returned after this function_
// →call.
// Specify the module using ModNum and the channel on the module using ChanNum._
// →Note that the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanBaselines (
    double *Baselines,                // returned baselines values
    double *TimeStamps,               // time stamp for each baseline value
    unsigned short NumBases,          // number of baseline values to read
    unsigned short ModNum,            // module number
    unsigned short ChanNum );         // channel number
```

```
// Ramp Offset DACs of a module and record the baselines
```

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```
// Use this function to execute the RAMP_OFFSETDACS control task run. Each Offset_
→DAC has 65536 steps, and the RAMP_OFFSETDACS control task ramps the DAC from 0_
→to 65335 with a step size of 64, i.e., a total of 1024 steps. At each DAC step,_
→the control task computes the baseline value as the representation of the signal_
→baseline and stores it in the DSP memory. After the control task is finished,_
→the stored baseline values are read out to the host computer and saved to a_
→binary file called "rampdacs.bin" in the form of IEEE 32-bit floating point_
→numbers. Users can then plot the baseline values vs. DAC steps to determine the_
→appropriate DAC value to be set in the DSP in order to bring the input signals_
→into the voltage range of the ADCs. However, this function is no longer needed_
→due to the introduction of function Pixie16AdjustOffsets.
// If ModNum is set to less than the total number of modules in the system, only_
→the module specified by ModNum will start the RAMP_OFFSETDACS control task run._
→But if ModNum is equal to the total number of modules in the system, e.g. there_
→are 5 modules in the chassis and ModNum = 5, then all modules in the system will_
→start the RAMP_OFFSETDACS control task run. Note that the modules are counted_
→starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RampOffsetDACS (
    double *DCValues,           // returned DC offset values
    unsigned short NumDCVals,   // number of DC values to read
    unsigned short ModNum );
```

```
// Execute special control tasks
// Use this function to call special control tasks. This may include programming_
→the Fippi or setting the DACs after downloading DSP parameters.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ControlTaskRun (
    unsigned short ModNum,       // Pixie module number
    unsigned short ControlTask, // Control task number
    unsigned int Max_Poll );   // Timeout control in unit of ms for control_
→task run
```

```
// Find the Baseline Cut values of a module
// Use this function to find the Baseline Cut value for one channel of a Pixie-16_
→module. The baseline cut value is then downloaded to the DSP, where baselines_
→are captured and averaged over time. The cut value would prevent a bad baseline_
→value from being used in the averaging process, i.e., if a baseline value is_
→outside the baseline cut range, it will not be used for computing the baseline_
→average. Averaging baselines over time improves energy resolution measurement.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BLcutFinder (
    unsigned short ModNum,       // Pixie module number
    unsigned short ChanNum,      // Pixie channel number
    unsigned int *BLcut );       // BLcut return value
```

```
// Find the exponential decay time of a channel
// Use this function to find the exponential decay time constant (Tau value) of_
→the detector or preamplifier signal that is connected to one channel of a Pixie-
→16 module. The found Tau value is returned via pointer *Tau.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16TauFinder (
    unsigned short ModNum,       // Pixie module number
    double *Tau );              // 16 returned Tau values, in [F]
```

```
// Write a MODULE level parameter to a module
// Use this function to write a module parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglModPar (
    char *ModParName,           // the name of the module parameter
    unsigned int ModParData,     // the module parameter value to be written to_
→the module
    unsigned short ModNum );    // module number
```

```
// Read a MODULE level parameter from a module
// Use this function to read a module parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglModPar (
    char *ModParName,           // the name of the module parameter
    unsigned int    *ModParData, // the module parameter value to be read from
    →the module
    unsigned short ModNum );   // module number
```

```
// Write a CHANNEL level parameter to a module
// Use this function to write a channel parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglChanPar (
    char *ChanParName,          // the name of the channel parameter
    double ChanParData,         // the channel parameter value to be written
    →to the module
    unsigned short ModNum,      // module number
    unsigned short ChanNum );   // channel number
```

```
// Read a CHANNEL level parameter from a module
// Use this function to read a channel parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanPar (
    char *ChanParName,          // the name of the channel parameter
    double *ChanParData,         // the channel parameter value to be read from
    →the module
    unsigned short ModNum,      // module number
    unsigned short ChanNum );   // channel number
```

```
// Read histogram data from a module
// Use this function to read out the histogram data from a Pixie-16 module's
// histogram memory. Before calling this function, the host code should allocate
// appropriate amount of memory to store the histogram data. The default histogram
// length is 32768. Histogram data are 32-bit unsigned integers.
// Specify the module using ModNum and the channel on the module using ChanNum.
// Note that both the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromModule (
    unsigned int    *Histogram,    // histogram data
    unsigned int    NumWords,     // number of words to be read out
    unsigned short ModNum,       // module number
    unsigned short ChanNum );    // channel number
```

```
// Read run statistics data from a module
// Use this function to read out statistics data from a Pixie-16 module. Before
// calling this function, the host code should allocate appropriate amount of
// memory to store the statistics data. The number of statistics data for each
// module is fixed at 448. Statistics data are 32-bit unsigned integers.
// Specify the module using ModNum. Note that the modules are counted starting at
// 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadStatisticsFromModule (
    unsigned int    *Statistics,   // run statistics data
    unsigned short ModNum );     // module number
```

```
// Read histogram data from a module and save to a file
// Use this function to read histogram data from a Pixie-16 module and save the
// data to a file. New data will be appended to the end of the file. So the same
// file name can be used for multiple modules and the data from each module will be
// stored in the order that this function is called.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveHistogramToFile (
    char *FileName,             // histogram data file name
    unsigned short ModNum);    // module number
```

```
// Parse a list mode data file to get events information
// Use this function to parse the list mode events in the list mode data file. The
// number of events for each module will be reported.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetModuleEvents (
    char *FileName, // the list mode data file name (with complete
    // path)
    unsigned int *ModuleEvents ); // receives number of events for each module
```

```
// Get detailed events information from a data file
// Use this function to retrieve the detailed information of each event in the
// list mode data file for the designated module. Before calling this function to
// get the individual events information, another function Pixie16GetModuleEvents
// should be called first to determine the number of events that have been recorded
// for each module. If the number of events for a given module is nEvents, a memory
// block *EventInformation should be allocated with a length of (nEvents*68):
// EventInformation = (unsigned long *)malloc(sizeof(unsigned long) * nEvents *
// 68);
// where 68 is the length of the information records of each event (energy,
// timestamps, etc.) and has the following structure.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetEventsInfo (
    char *FileName, // the list mode data file name (with
    // complete path)
    unsigned int *EventInformation, // to hold event information
    unsigned short ModuleNumber); // the module whose events are to be
    // retrieved
```

```
// Read trace data from a list mode data file
// Use this function to retrieve list mode trace from a list mode data file. It
// uses the trace length and file location information obtained from function
// Pixie16GetEventsInfo for the selected event.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadListModeTrace (
    char *FileName, // list mode data file name
    unsigned short *Trace_Data, // list mode trace data (16-bit words)
    unsigned short NumWords, // number of 16-bit words to be read out
    unsigned int FileLocation); // the location of the trace in the file
```

```
// Read histogram data from a histogram data file
// Use this function to read histogram data from a histogram data file. Before
// calling this function, the host code should allocate appropriate amount of
// memory to store the histogram data. The default histogram length is 32768.
// Histogram data are 32-bit unsigned integers.
// Specify the module using ModNum and the channel on the module using ChanNum.
// Note that both the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromFile (
    char *FileName, // the histogram data file name (with complete
    // path)
    unsigned int *Histogram, // histogram data
    unsigned int NumWords, // number of words to be read out
    unsigned short ModNum, // module number
    unsigned short ChanNum); // channel number
```

```
// Read DSP parameters from modules and save to a file
// Use this function to save DSP parameters to a settings file. It will first read
// the values of DSP parameters on each Pixie-16 module and then write them to the
// settings file. Each module has exactly 1280 DSP parameter values (32-bit
// unsigned integers), and depending on the value of PRESET_MAX_MODULES (defined in
// pixie16app_defs.h), the settings file should have exactly (1280 * PRESET_MAX_
// MODULES * 4) bytes when stored on the computer hard drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveDSPParametersToFile (
    char *FileName ); // the DSP parameters file name (with complete
    // path)
```

```
// Load DSP parameters to modules from a file
// Use this function to read DSP parameters from a settings file and then download_
// the settings to Pixie-16 modules that are installed in the system. Each module_
// has exactly 1280 DSP parameter values (32-bit unsigned integers), and depending_
// on the value of PRESET_MAX_MODULES (defined in pixie16app_defs.h), the settings_
// file should have exactly (1280 * PRESET_MAX_MODULES * 4) bytes when stored on_
// the computer hard drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16LoadDSPParametersFromFile (
    char *FileName ); // the DSP parameters file name (with complete_
// path)
```

```
// Copy DSP parameters from a module to others
// Use this function to copy DSP parameters from one module to the others that are_
// installed in the system.
// BitMask is bit pattern which designates which items should be copied from the_
// source module to the destination module(s).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CopyDSPParameters (
    unsigned short BitMask, // copy items bit mask
    unsigned short SourceModule, // source module
    unsigned short SourceChannel, // source channel
    unsigned short *DestinationMask ); // the destination module and channel_
// bit mask
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadMSGFile (
    char *ReturnMsgStr );
```

```
// Convert a decimal into IEEE 32-bit floating point number
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API Decimal2IEEEFloating(double_
//DecimalNumber);

// Convert an IEEE 32-bit floating point number to a decimal
PIXIE16APP_EXPORT double PIXIE16APP_API IEEEFloating2Decimal(unsigned int_
//IEEEfloatingNumber);
```

```
// Read data from external FIFO and save to a file
// Use this function to read data from the external FIFO of a module. This_
// function can only be used for Pixie-16 Revision-B, C, and D modules.
// This function first checks the status of the external FIFO of a Pixie-16 module,
// and if there are data in the external FIFO, this function then reads list mode_
// data (32-bit unsigned integers) from the external FIFO. So this function_
// essentially encapsulates both functions Pixie16CheckExternalFIFOStatus and_
// Pixie16ReadDataFromExternalFIFO within one function. The number of words that_
// are read from the external FIFO is recorded in variable*nFIFOWords.
// The function also expects setting the value of a variable called "EndOfRunRead"
// to indicate whether this read is at the end of a run (1) or during the run (0)._
// This is necessary since the external FIFO needs special treatment when the host_
// reads the last few words from the external FIFO due to its pipelined structure.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveExternalFIFODataToFile (
    char *FileName, // list mode data file name
    unsigned int *nFIFOWords, // number of words read from external FIFO
    unsigned short ModNum, // module number
    unsigned short EndOfRunRead); // indicator whether this is the end of run_
//read
```

```
// Read from or write to registers on a module
// Use this function to read data from or write data to a register in a Pixie-16_
// module.
// Specify the module using ModNum. Note that the modules are counted starting at_
// 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RegisterIO (
```

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```

unsigned short ModNum,           // the Pixie module to communicate to
unsigned int address,          // register address
unsigned short direction,       // either MOD_READ or MOD_WRITE
unsigned int *value );        // holds or receives the data

// Read Control & Status Register value from a module
// Use this function to read the host Control & Status Register (CSR) value.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadCSR (
    unsigned short ModNum,
    unsigned int *CSR );

// Write to Control & Status Register in a module
// Use this function to write a value to the host Control & Status Register (CSR).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteCSR (
    unsigned short ModNum,
    unsigned int CSR );

// Check status of external FIFO of a module
// Use this function to check the status of the external FIFO of a Pixie-16 module.
// While a list mode data acquisition run is in progress. The function returns the
// number of words (32-bit) that the external FIFO currently has. If the number of
// words is greater than a user-set threshold, function
// Pixie16ReadDataFromExternalFIFO can then be used to read the data from the
// external FIFO. The threshold can be set by the user to either minimize reading
// overhead or to read data out of the FIFO as quickly as possible.
// *nFIFOWords returns the number of 32-bit words that the external FIFO currently
// has.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckExternalFIFOStatus (
    unsigned int *nFIFOWords,
    unsigned short ModNum );

// Read data from external FIFO of a module
// Use this function to read data from the external FIFO of a module. This
// function can only be used for Pixie-16 Revision-B, C, and D modules.
// This function reads list mode data from the external FIFO of a Pixie-16 module.
// The data are 32-bit unsigned integers. Normally, function
// Pixie16CheckExternalFIFOStatus is called first to see how many words the
// external FIFO currently has, then this function is called to read the data from
// the FIFO.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadDataFromExternalFIFO (
    unsigned int *ExtFIFO_Data, // To receive the external FIFO data
    unsigned int nFIFOWords,   // number of words to read from external FIFO
    unsigned short ModNum ); // module number

PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeFastFiltersOffline (
    char             *FileName,           // the list mode data file name (with
    // complete path)
    unsigned short ModuleNumber,        // the module whose events are to be
    // analyzed
    unsigned short ChannelNumber,       // the channel whose events are to be
    // analyzed
    unsigned int   FileLocation,         // the location of the trace in the file
    unsigned short RcdTraceLength,      // recorded trace length
    unsigned short *RcdTrace,           // recorded trace
    double          *fastfilter,        // fast filter response
    double          *cfd );            // cfd response

```

```

PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with_
→complete path)
    unsigned short ModuleNumber,      // the module whose events are to be_
→analyzed
    unsigned short ChannelNumber,     // the channel whose events are to be_
→analyzed
    unsigned int   FileLocation,      // the location of the trace in the file
    unsigned short RcdTraceLength,    // recorded trace length
    unsigned short *RcdTrace,         // recorded trace
    double         *slowfilter );    // slow filter response

```

```

// Add by Hongyi Wu
PIXIE16APP_EXPORT int PIXIE16APP_API HongyiWuPixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with_
→complete path)
    unsigned short ModuleNumber,      // the module whose events are to be_
→analyzed
    unsigned short ChannelNumber,     // the channel whose events are to be_
→analyzed
    unsigned int   FileLocation,      // the location of the trace in the file
    unsigned short RcdTraceLength,    // recorded trace length
    unsigned short *RcdTrace,         // recorded trace
    double         *slowfilter,        // slow filter response
    unsigned int   bl,
    double         sl,
    double         sg,
    double         tau,
    int            sfr,
    int            pointtbl );

```

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